



WINNER 2025
Bi
BEST OF
INDUSTRY
A W A R D

HOT MOUNTING PRESS

QPRESS 40

qatm.com/products/mounting/qpress-40/

QPRESS 40

HOT MOUNTING PRESS

The Qpress 40 is a robust hot mounting press for exceptionally fast mounting processes. Thanks to its powerful heating/cooling technology and speed-optimized temperature control, the mounting time is significantly reduced.

The robust bayonet closure, optional simple one-button operation, and clear parameter display on the 7" industrial touchscreen ensure convenient and intuitive use. The connectable, patented extraction system sets new standards in occupational safety for the user.

The new Qpress 40 hot mounting press has been developed to meet the increasing demands for precision, efficiency and user-friendliness. The Qpress 40 is characterized by innovative features that not only optimize work processes, but also improve the quality of the results.




QPRESS 40

PRODUCT ADVANTAGES

- Short mounting times (up to 35% time savings)
- 7" industrial touchscreen with innovative control software and optional one-button operation
- Compact and high-quality design with powder-coated metal housing
- Robust bayonet quick-release closure
- Dust Guard: extraction of granule dust during filling (patented)
- High sample throughput possible due to double mounting

QPRESS 40

FEATURES

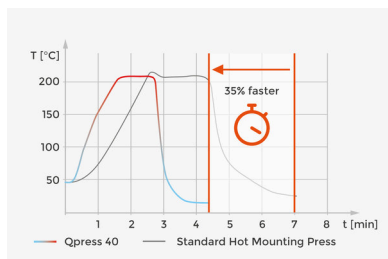**SPEED-OPTIMIZED HEATING AND COOLING TECHNOLOGY**

Shorter mounting times are achieved by optimized heat transfer between heating elements and mould, fast temperature control to precisely reach the target temperature, and a preheating function to shorten the heat-up time. The new technology can reduce the mounting time by up to 35%.

**SIMPLE ONE-BUTTON OPERATION**

The highest possible operating safety is achieved with the butler function. In just a few steps, the entire mounting process can be controlled by pressing the multifunction button. PLUS software has a method database with pre-installed recipes and stores more than 128 user-defined mounting methods.

**EXTRACTION OF FINE DUST**

The patented and optional integrable suction device (Dust Guard) protects the user from inhaling dust that can escape when the mounting resin is filled into the mould. In the PLUS software version, extraction is automatically controlled by the software.

**EASY TO CLEAN**

The suction nozzle attached to the back of the machine is removable, making cleaning very convenient. The mould and the two plungers can be removed for cleaning in just a few simple steps.

QPRESS 40

FEATURES CONTINUED



ROBUST BAYONET CLOSURE

The robust and ergonomic bayonet closure fits comfortably in the hand and allows quick and easy opening and closing of the press cylinder.



TALL PRESS CYLINDER FOR DOUBLE MOUNTING

The tall press cylinder allows two samples to be mounted on top of each other simultaneously. This saves valuable mounting time and doubles sample throughput. Additionally, this mounting method is particularly time- and energy-efficient.




GRANULATE COMPARTMENT

The machine includes a granule compartment for easy pick-up and disposal of excess mounting granules.

QPRESS 40

FUNCTIONAL PRINCIPLE

The new design of the Qpress 40 ensures the shortest possible heating and cooling times. The heating elements quickly melt the mounting resin thanks to optimized heat transfer and temperature control in the mould. During this process, the granules are compressed to the preset mould pressure. Once the set temperature is reached, the press mould is cooled as quickly as possible after a short holding time, allowing the hardened sample to be removed for further preparation. This not only saves time but also reduces energy requirements.

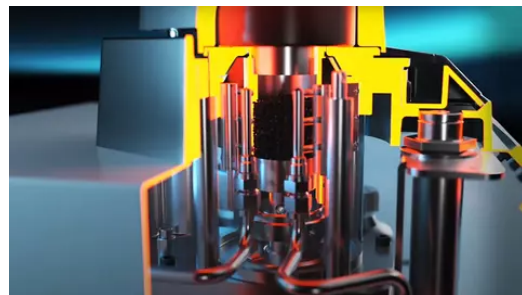

QPRESS 40

MOUNTING IN A FEW STEPS

**1 | PLACE SAMPLE**

Place sample.

**2 | FILL IN MOUNTING RESIN**

Fill in mounting resin.

**3 | LOCK THE CLOSURE**

Lock the bayonet quick-release closure.

**4 | START PROCESS**

Start process.

**5 | REMOVE THE SAMPLE**

Open the closure and remove the sample.

QPRESS 40


SOFTWARE

The Qpress 40 software is available in two variants. The standard Go software allows for quick setting of process parameters and starting of the mounting process. The Plus software variant offers additional features such as a method database, butler function, and automatic control of extraction. The Plus software also includes the 2-step mode for sensitive samples. You can upgrade to Advanced Software Plus at any time.

FUNCTION	GO (Standard)	PLUS (Advanced)
Setting of the mounting parameters via industrial touchscreen (heating/ cooling times, temperature, pressure)	✓	✓
Speed-optimized temperature control	✓	✓
Possibility of double mounting	✓	✓
User account control	✓	✓
Status bar	✓	✓
Adjustable cooling capacity	✓	✓
Automatic descaling function	✓	✓
Butler function (simple one-button operation)	✗	✓
Dust Guard: Extraction of fine dust (controlled via software)	✗	✓
Method database (saving and loading programs)	✗	✓
2-Step Mode (for sensitive samples)	✗	✓


QPRESS 40

SOFTWARE IN PRACTICE

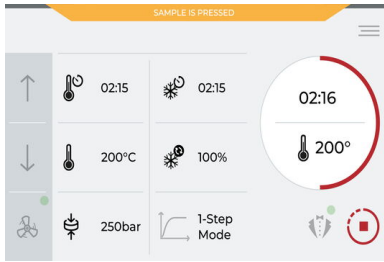

CLEAR SOFTWARE

The intuitive software displays all important mounting parameters on the robust 7" industrial touchscreen. Set heating and cooling times, cooling capacity, temperature and press mould pressure. The software displays the press status; User Account Control allows settings to be write-protected.

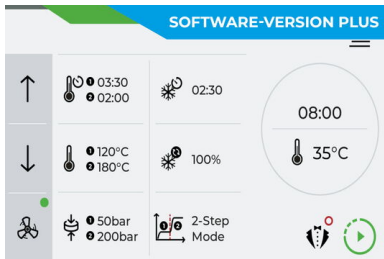

2-STEP PRINT MODE

For sensitive samples (sheets, thin-walled pipes, coated materials), a 2-step mode permits independent control of temperature and pressure. The mounting resin melts at low pressure first; full pressure is then applied for compaction.



SOFTWARE-SUPPORTED DESCALING FUNCTION

The machine features an automatic descaling function. Simply insert a descaling tablet and select the software function to start the descaling process.

QPRESS 40

BEST OF INDUSTRY AWARDS 2025

The Qpress 40 hot mounting press was awarded the Best of Industry Award 2025 in the category "Best of Laboratory and Analysis Technology". The award recognizes innovative technologies that make a lasting improvement to everyday lab work and offer measurable added value for research, diagnostics, and applications.



QPRESS 40

ACCESSORIES



CIRCULATING COOLING

- | Powerful cooling
- | Resource-saving circulation system instead of fresh water
- | Low maintenance due to closed-circuit coolant circulation, reducing limescale deposits
- | Can connect up to two hot mounting presses



MOULDS

Easily interchangeable moulds are available with or without chamfer in these sizes: Ø 25 mm, Ø 30 mm, Ø 40 mm, Ø 1", Ø 1 1/4", Ø 1 1/2".



SUCTION

The optional extraction device allows direct dust extraction when filling the mounting resin and easy cleaning with a removable suction hose. It connects to up to two mounting machines and is controlled via Qpress 40 software. It includes a replaceable vacuum cleaner bag and an additional fine dust filter.

QPREP

PERFECT MOUNTING – CONSUMABLES

Mounting is more reliable, safer and more efficient when performed with the right equipment: Knowledge, exercising care and the right choice of consumables and equipment are crucial.

TO THE PRODUCT RANGE →


qatm.com/products/consumables/mounting/hot-mounting/



QPRESS 40

TECHNICAL DATA

Mould Ø (metric)	Ø 25 mm, Ø 30 mm, Ø 40 mm
Mould Ø (inch)	Ø 1", Ø 1 1/4", Ø 1 1/2"
Max. pressure in mould cylinder Ø 40 mm / Ø 1.5"	250 bar
Max. pressure in mould cylinder < Ø 40 mm	350 bar
Adjustable temperature	20 - 200°C
Drive power (hydraulics)	0.36 kW
Heating capacity	1.6 kW
Connected load	2.0 kW
Heating time (adjustable)	0 - 30 min; in 15 s steps
Cooling time (adjustable)	0 - 30 min; in 15 s steps
Dimensions (WxHxD)	~ 365 x 433 x 538mm
Weight	~ 41 kg

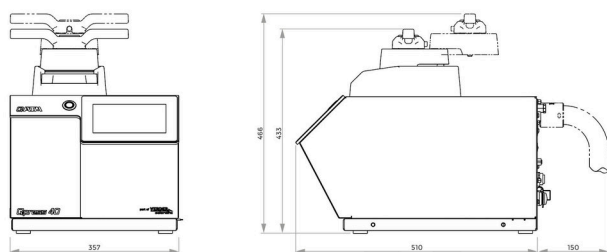

QPRESS 40

DOWNLOADS

PRODUCT DATA SHEET → qatm.com/files/450085/product-data-sheet-qpress-40.pdf

PRODUCT BROCHURE → qatm.com/files/446721/qpress-40-qpress-50-qmount.pdf

QPRESS 40

ORDERING INFORMATION

BASIC EQUIPMENT

MACHINE (1 item required)

M0770020	Qpress 40 basic unit, 200-240 V 50/60 Hz (1Ph/N/PE)
M0770021	Qpress 40 basic unit, 100-120 V 50/60 Hz (1Ph/N/PE)

SOFTWARE (can be ordered separately)

Z0770220	Qpress 40 upgrade for software version PLUS
-----------------	---

MOULD ASSEMBLY (can be ordered separately)

Z0770200	Qpress 40 mould assembly Ø 25 mm
Z0770201	Qpress 40 mould assembly Ø 1"
Z0770202	Qpress 40 mould assembly Ø 30 mm
Z0770203	Qpress 40 mould assembly Ø 1 1/4"
Z0770204	Qpress 40 mould assembly Ø 1 1/2"
Z0770205	Qpress 40 mould assembly Ø 40 mm
Z0770210	Qpress 40 mould assembly Ø 25 mm with chamfer
Z0770211	Qpress 40 mould assembly Ø 1" with chamfer
Z0770212	Qpress 40 mould assembly Ø 30 mm with chamfer
Z0770213	Qpress 40 mould assembly Ø 1 1/4" with chamfer
Z0770214	Qpress 40 mould assembly Ø 1 1/2" with chamfer
Z0770215	Qpress 40 mould assembly Ø 40 mm with chamfer

SPACER (can be ordered separately)

00770306	Spacer Ø 25 mm for double mounts
00770307	Spacer Ø 30 mm for double mounts
00770308	Spacer Ø 40 mm for double mounts
00770309	Spacer Ø 1" for double mounts
00770310	Spacer Ø 1 1/4" for double mounts
00770311	Spacer Ø 1 1/2" for double mounts
00770312	Spacer Ø 25 mm with chamfer for double mounts
00770313	Spacer Ø 30 mm with chamfer for double mounts
00770314	Spacer Ø 40 mm with chamfer for double mounts
00770315	Spacer Ø 1" with chamfer for double mounts
00770316	Spacer Ø 1 1/4" with chamfer for double mounts
00770317	Spacer Ø 1 1/2" with chamfer for double mounts

OPTIONS (can be ordered separately)

95014560	Starter Kit Hot Mounting
Z0800001	Connection set
Z0770226	Cooling unit, 200-240V 50/60Hz (1Ph/N/PE)
Z0770227	Cooling unit, 100-120V 50/60Hz (1Ph/N/PE)
Z0770225	Connection set vacuum cleaner second machine
Z0770223	Vacuum cleaner, 200-240V 50/60Hz (1Ph/N/PE)
Z0770224	Vacuum cleaner, 100-120V 50/60Hz (1Ph/N/PE)

QPRESS 40

CONTACT US →

qatm.com/contact/contact-form/?prefilltext=contact_form.request&product=A01B07A

REQUEST A QUOTE →

qatm.com/products/mounting/qpress-40/order-data/



WINNER 2025
BI
BEST OF
INDUSTRY
AWARD



qatm.com/products/mounting/qpress-40/