



WINNER 2025
Bi
BEST OF
INDUSTRY
A W A R D

ПРЕССЫ ГОРЯЧЕГО ПРЕССОВАНИЯ

QPRESS 40

qatm.com/ru/products/mounting/qpress-40/

QPRESS 40

ПРЕССЫ ГОРЯЧЕГО ПРЕССОВАНИЯ

The Qpress 40 is a robust hot mounting press for exceptionally fast mounting processes. Thanks to its powerful heating/cooling technology and speed-optimized temperature control, the mounting time is significantly reduced.

The robust bayonet closure, optional simple one-button operation, and clear parameter display on the 7" industrial touchscreen ensure convenient and intuitive use. The connectable, patented extraction system sets new standards in occupational safety for the user.

The new Qpress 40 hot mounting press has been developed to meet the increasing demands for precision, efficiency and user-friendliness. The Qpress 40 is characterized by innovative features that not only optimize work processes, but also improve the quality of the results.




QPRESS 40

ПРЕИМУЩЕСТВА ОБОРУДОВАНИЯ

- | Short mounting times (up to 35% time savings)
- | 7" industrial touchscreen with innovative control software and optional one-button operation
- | Compact and high-quality design with powder-coated metal housing
- | Robust bayonet quick-release closure
- | Dust Guard: extraction of granule dust during filling (patented)
- | High sample throughput possible due to double mounting

QPRESS 40

ОСОБЕННОСТИ

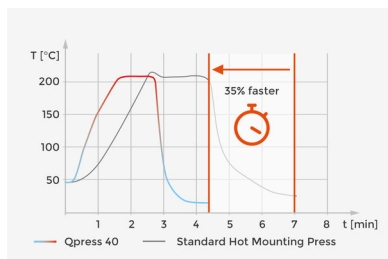

SPEED-OPTIMIZED HEATING AND COOLING TECHNOLOGY

Shorter mounting times are achieved by optimized heat transfer between heating elements and mould, fast temperature control to precisely reach the target temperature, and a preheating function to shorten the heat-up time. The new technology can reduce the mounting time by up to 35%.



SIMPLE ONE-BUTTON OPERATION

The highest possible operating safety is achieved with the butler function. In just a few steps, the entire mounting process can be controlled by pressing the multifunction button. PLUS software has a method database with pre-installed recipes and stores more than 128 user-defined mounting methods.



EXTRACTION OF FINE DUST

The patented and optional integrable suction device (Dust Guard) protects the user from inhaling dust that can escape when the mounting resin is filled into the mould. In the PLUS software version, extraction is automatically controlled by the software.



EASY TO CLEAN

The suction nozzle attached to the back of the machine is removable, making cleaning very convenient. The mould and the two plungers can be removed for cleaning in just a few simple steps.

QPRESS 40

ОСОБЕННОСТИ: ПРОДОЛЖЕНИЕ



ROBUST BAYONET CLOSURE

The robust and ergonomic bayonet closure fits comfortably in the hand and allows quick and easy opening and closing of the press cylinder.



TALL PRESS CYLINDER FOR DOUBLE MOUNTING

The tall press cylinder allows two samples to be mounted on top of each other simultaneously. This saves valuable mounting time and doubles sample throughput. Additionally, this mounting method is particularly time- and energy-efficient.




GRANULATE COMPARTMENT

The machine includes a granule compartment for easy pick-up and disposal of excess mounting granules.

QPRESS 40

ПРИНЦИП РАБОТЫ

The new design of the Qpress 40 ensures the shortest possible heating and cooling times. The heating elements quickly melt the mounting resin thanks to optimized heat transfer and temperature control in the mould. During this process, the granules are compressed to the preset mould pressure. Once the set temperature is reached, the press mould is cooled as quickly as possible after a short holding time, allowing the hardened sample to be removed for further preparation. This not only saves time but also reduces energy requirements.

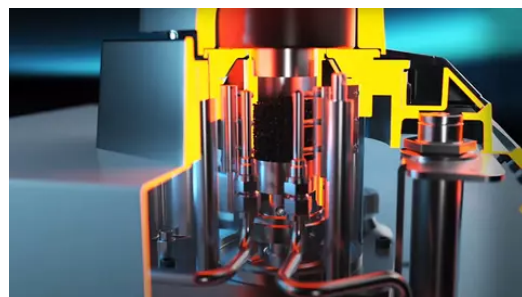

QPRESS 40

MOUNTING IN A FEW STEPS

**1 | PLACE SAMPLE**

Place sample.

**2 | FILL IN MOUNTING RESIN**

Fill in mounting resin.

**3 | LOCK THE CLOSURE**

Lock the bayonet quick-release closure.

**4 | START PROCESS**

Start process.

**5 | REMOVE THE SAMPLE**

Open the closure and remove the sample.

QPRESS 40


ПРОГРАММНОЕ ОБЕСПЕЧЕНИЕ

The Qpress 40 software is available in two variants. The standard Go software allows for quick setting of process parameters and starting of the mounting process. The Plus software variant offers additional features such as a method database, butler function, and automatic control of extraction. The Plus software also includes the 2-step mode for sensitive samples. You can upgrade to Advanced Software Plus at any time.

FUNCTION	GO (Standard)	PLUS (Advanced)
Setting of the mounting parameters via industrial touchscreen (heating/cooling times, temperature, pressure)	✓	✓
Speed-optimized temperature control	✓	✓
Possibility of double mounting	✓	✓
User account control	✓	✓
Status bar	✓	✓
Adjustable cooling capacity	✓	✓
Automatic descaling function	✓	✓
Butler function (simple one-button operation)	✗	✓
Dust Guard: Extraction of fine dust (controlled via software)	✗	✓
Method database (saving and loading programs)	✗	✓
2-Step Mode (for sensitive samples)	✗	✓


QPRESS 40

ПРОГРАММНОЕ ОБЕСПЕЧЕНИЕ: ПРАКТИКА

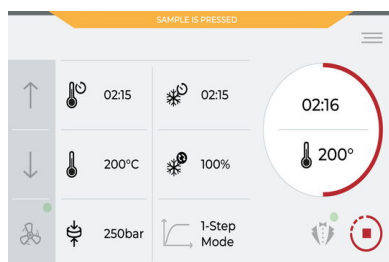

CLEAR SOFTWARE

The intuitive software displays all important mounting parameters on the robust 7" industrial touchscreen. Set heating and cooling times, cooling capacity, temperature and press mould pressure. The software displays the press status; User Account Control allows settings to be write-protected.

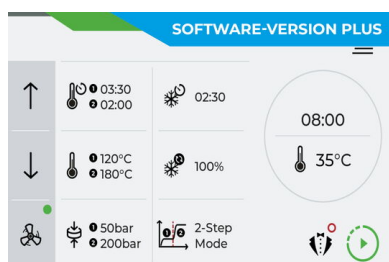

2-STEP PRINT MODE

For sensitive samples (sheets, thin-walled pipes, coated materials), a 2-step mode permits independent control of temperature and pressure. The mounting resin melts at low pressure first; full pressure is then applied for compaction.



SOFTWARE-SUPPORTED DESCALING FUNCTION

The machine features an automatic descaling function. Simply insert a descaling tablet and select the software function to start the descaling process.

QPRESS 40

BEST OF INDUSTRY AWARDS 2025

The Qpress 40 hot mounting press was awarded the Best of Industry Award 2025 in the category "Best of Laboratory and Analysis Technology". The award recognizes innovative technologies that make a lasting improvement to everyday lab work and offer measurable added value for research, diagnostics, and applications.



QPRESS 40

АКСЕССУАРЫ



CIRCULATING COOLING

- | Powerful cooling
- | Resource-saving circulation system instead of fresh water
- | Low maintenance due to closed-circuit coolant circulation, reducing limescale deposits
- | Can connect up to two hot mounting presses



MOULDS

Easily interchangeable moulds are available with or without chamfer in these sizes: Ø 25 mm, Ø 30 mm, Ø 40 mm, Ø 1", Ø 1 1/4", Ø 1 1/2".



SUCTION

The optional extraction device allows direct dust extraction when filling the mounting resin and easy cleaning with a removable suction hose. It connects to up to two mounting machines and is controlled via Qpress 40 software. It includes a replaceable vacuum cleaner bag and an additional fine dust filter.

QPREP

PERFECT MOUNTING – CONSUMABLES BY QPREP

запрессовка надежнее, безопаснее и эффективнее, если она выполняется с использованием правильного оборудования: знания, осторожность и правильный выбор расходных материалов и оборудования имеют решающее значение.

К КАТАЛОГУ →


qatm.com/ru/products/consumables/mounting/hot-mounting/



QPRESS 40

ХАРАКТЕРИСТИКИ

Плесень $\varnothing$ (метрика)	$\varnothing$ 25 мм, $\varnothing$ 30 мм, $\varnothing$ 40 мм
Mould $\varnothing$ (inch)	$\varnothing$ 1", $\varnothing$ 1 1/4", $\varnothing$ 1 1/2"
Макс. давление в цилиндре формы $\varnothing$ 40 мм / $\varnothing$ 1,5"	250 Бар
Макс. давление в цилиндре формы < $\varnothing$ 40 мм	350 бар
Регулируемая температура	20 - 200°C
Мощность привода (гидравлика)	0.36 кВт
Мощность отопления	1,6 кВт
Подключенная нагрузка	2,0 кВт
Время нагрева (регулируемое)	0 - 30 мин; в 15 с шагов
Время охлаждения (регулируемое)	0 - 30 мин; в 15 с шагов
Размер (ВхШхГ)	~ 365 x 433 x 538 мм
Вес	~ 41 кг

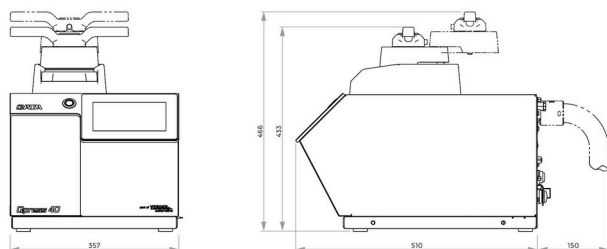

QPRESS 40

СКАЧАТЬ

ПАСПОРТ ИЗДЕЛИЯ QPRESS 40 → qatm.com/files/450536/pasport-izdelia-qpress-40.pdf

QPRESS 40 / QPRESS 50 / QMOUNT → qatm.com/files/446721/qpress-40-qpress-50-qmount.pdf

QPRESS 40

ИНФОРМАЦИЯ ДЛЯ ЗАКАЗА

БАЗОВОЕ ОБОРУДОВАНИЕ

МАШИНА (необходима - 1 шт.)

M0770020	Qpress 40 basic unit, 200-240 V 50/60 Hz (1Ph/N/PE)
----------	---

M0770021	Qpress 40 basic unit, 100-120 V 50/60 Hz (1Ph/N/PE)
----------	---

SOFTWARE (can be ordered separately)

Z0770220	Qpress 40 upgrade for software version PLUS
----------	---

ПРЕСС-ФОРМА (может заказываться отдельно)

Z0770200	Qpress 40 mould assembly Ø 25 mm
----------	----------------------------------

Z0770201	Qpress 40 mould assembly Ø 1"
----------	-------------------------------

Z0770202	Qpress 40 mould assembly Ø 30 mm
----------	----------------------------------

Z0770203	Qpress 40 mould assembly Ø 1 1/4"
----------	-----------------------------------

Z0770204	Qpress 40 mould assembly Ø 1 1/2"
----------	-----------------------------------

Z0770205	Qpress 40 mould assembly Ø 40 mm
----------	----------------------------------

Z0770210	Qpress 40 mould assembly Ø 25 mm with chamfer
----------	---

Z0770211	Qpress 40 mould assembly Ø 1" with chamfer
----------	--

Z0770212	Qpress 40 mould assembly Ø 30 mm with chamfer
----------	---

Z0770213	Qpress 40 mould assembly Ø 1 1/4" with chamfer
----------	--

Z0770214	Qpress 40 mould assembly Ø 1 1/2" with chamfer
----------	--

Z0770215	Qpress 40 mould assembly Ø 40 mm with chamfer
----------	---

ДИСТАНЦИОННОЕ КОЛЬЦО (может заказываться отдельно)

00770306	Spacer Ø 25 mm for double mounts
----------	----------------------------------

00770307	Spacer Ø 30 mm for double mounts
----------	----------------------------------

00770308	Spacer Ø 40 mm for double mounts
----------	----------------------------------

00770309	Spacer Ø 1" for double mounts
----------	-------------------------------

00770310	Spacer Ø 1 1/4" for double mounts
----------	-----------------------------------

00770311	Spacer Ø 1 1/2" for double mounts
----------	-----------------------------------

00770312	Spacer Ø 25 mm with chamfer for double mounts
----------	---

00770313	Spacer Ø 30 mm with chamfer for double mounts
----------	---

00770314	Spacer Ø 40 mm with chamfer for double mounts
----------	---

00770315	Spacer Ø 1" with chamfer for double mounts
----------	--

00770316	Spacer Ø 1 1/4" with chamfer for double mounts
----------	--

00770317	Spacer Ø 1 1/2" with chamfer for double mounts
----------	--

ОПЦИИ (могут заказываться отдельно)

95014560	Starter Kit Hot Mounting
----------	--------------------------

Z0800001	Комплект для подключения
----------	--------------------------

Z0770226	Cooling unit, 200-240V 50/60Hz (1Ph/N/PE)
----------	---

Z0770227	Cooling unit, 100-120V 50/60Hz (1Ph/N/PE)
----------	---

Z0770225	Connection set vacuum cleaner second machine
----------	--

Z0770223	Vacuum cleaner, 200-240V 50/60Hz (1Ph/N/PE)
----------	---

Z0770224	Vacuum cleaner, 100-120V 50/60Hz (1Ph/N/PE)
----------	---

QPRESS 40

СВЯЖИТЕСЬ С НАМИ →

qatm.com/ru/contact/contact-form/?prefilltext=contact_form.request&product=A01B07A

ЗАПРОС ПРЕДЛОЖЕНИЯ →

qatm.com/ru/products/mounting/qpress-40/order-data/



WINNER 2025
BI
BEST OF
INDUSTRY
AWARD

QATM
a **VERDER** company

qatm.com/ru/products/mounting/qpress-40/