



PRECISIE DOORSLIJPMACHINE

QCUT 150 A

De Qcut 150 A is een flexibele doorslijpmachine die gebruikt wordt met slijpschijven tot diameter 203 mm/8".

The compact benchtop model performs automatic as well as manual cuts with the ergonomic lever. The optional sample holder arm with built-in force sensor ensures low-deformation cuts due to a precise force-controlled feed rate.

PRODUCTVOORDELEN

- | Tot 5 verschillende modi: manueel, automatische aanvoer, kracht-gecontroleerde aanvoer, slijpmodus, snijden met tafel bevestiging (afhankelijk van machine-configuratie)
- | Manuele dwarsaanvoer
- | Functie "zone-slijpen"
- | 4 verschillende koelsystemen beschikbaar
- | Slijpsteen-aftrekker
- | Intuïtieve QATM software met geïntegreerde databank



[Klik om video te bekijken](#)

DIT INSTRUMENT IS IDEAAL
VOOR GEOLOGISCHE EN
MINERALOGISCHE
TOEPASSINGEN

PRECISIE DOORSLIJPMACHINE QCUT 150 A

KENMERKEN




KRACHT-GESTUURD DOORSLIJPEN

- | Integrated force sensor in the sample holder arm precisely regulates the feed rate via the set cutting force
- | Constant cutting force throughout the process for clean cutting even of edge-hardened components and composite materials



TAFEL BEVESTIGING

- | Table attachment can be mounted without tools
- | Quick and manual sectioning of large samples, PCBs, or plate material

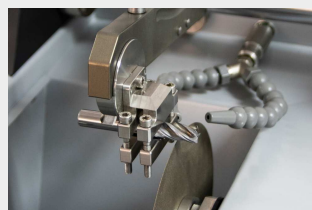

NIEUW: ZONE-SLIJPEN

- | Individual adjustment of feed rate or cutting force for single sections of the sample
- | Reduced material deformation and heating
- | Accelerated cutting process




ALL PARAMETERS AT A GLANCE

- | Robust 4.3" industrial touchscreen
- | Intuitive software with modern interface displaying all important cutting parameters
- | Extensive program storage and user account control



SLIJP-MODUS

- | Target preparation: In grinding mode, samples can be ground manually or semi-automatically using a cup wheel
- | Ideal for thin sections and microelectronics

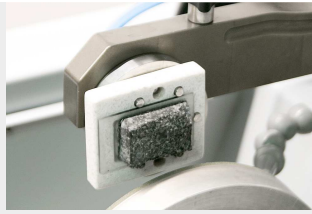

VACUUM CLAMPING UNIT

- | Glass slides for thin sections can be fixed easily and without slipping

REDUCES SUBSEQUENT GRINDING PROCESSES

The use of the adaptive rotation axis leads to an even more uniform cutting pattern and reduces the risk of the cut-off wheel runout. This reduces the effort required for later plane-parallel grinding of the sample.



PRECISIE DOORSLIJPMACHINE QCUT 150 A

SUPERIORITY IN DETAIL

1. Sample holder arm

- | Tool-free quick-clamping system for mounting :
- | Precision clamping devices
- | Rotation device
- | Vacuum clamping unit
- | Max. sample size: Ø

5. Transparent hood

- | Swivel hood with free access to the cutting chamber
- | Integrated safety switch

6. Z-axis

- | Manual cross feed

40 mm

2. Y-axis

- | Manual or automatic feed
- | Travel: 120 mm
- | Feed rate: 0.01–20 mm/s (automatic)
- | Cutting force variably adjustable: 0–10 N (optional)

3. Dressing device

- | For sharpening diamond cutting and grinding wheels
- | Can be operated manually from the outside

4. Robust machine body

- | Sturdy and corrosion-resistant stainless steel housing
- | Powder-coated, for easier cleaning

- | Travel: 25 mm
- | Position shown on the display
- | Positioning accuracy: 0.005 mm

7. Cut-off Wheel drive

- | Cut-off wheel size: Ø 75–203 mm / 8"
- | Drive power: 0.37 kW (SI)

8. Ergonomic Control panel

- | 4.3" industrial touchscreen
- | Intuitive interface displaying process data
- | Storage of up to 200 user-defined cutting methods
- | User account management





**UITGEBREIDE REEKS
KLEMHULPSTUKKEN**

With the help of a suitable clamping tool from the extensive range of accessories, mounted or non-mounted samples up to a size of Ø 40 mm can be clamped in the cut-off machine comfortably

OPTIMAL ACCESSORIES FOR QCUT 150 A AND 200 A

COOLING SYSTEMS

Optimal cooling is essential during sectioning. Excessive heat generation during the sectioning process damages the specimen and increases the wear of the cut-off wheel. The cooling and corrosion protection agent dissipates the heat and simultaneously removes the chips from the cut-off chamber.



15L RECIRCULATION COOLING SYSTEM

- | Portable system with filter cartridge and lid
- | Small footprint required
- | Easy to integrate into existing laboratory furniture
- | Two-chamber system for filtration



45L RECIRCULATION COOLING SYSTEM (MOBILE)

- | Mobile unit with inlet strainer, filter bag and overflow
- | Can be placed under a lab bench or floor-to-ceiling cabinet
- | Expandable with optional pump module for cut-off chamber flushing
- | Two-chamber system for filtration with removable separating slide



45L RECIRCULATION COOLING SYSTEM (SYSTEM LAB)

- | Mobile unit with retractable mechanism
- | Can be conveniently integrated into QATM system laboratory cabinets and cabinets with undershelves
- | Expandable with optional pump module for cut-off chamber flushing
- | Inlet strainer, filter bag and overflow
- | Two-chamber system for filtration with removable separating slide



CUT-OFF WHEEL IMMERSION COOLING

- | Equipment to seal and fill cutting chamber with coolant

- | Incl. drain hose and shut-off valve for easy emptying
- | Sample collector

PERFECT CUTTING - CUT-OFF WHEELS FROM QATM

Cutting is more reliable, safe and efficient when processed with the right ingredients: knowledge, care and the right choice of consumables and instruments. The Qprep range includes specially developed cut-off wheels for a wide variety of materials and requirements.

- | Low thermal influence on the specimen during the cutting process
- | Wet abrasive cutting achieves very fine surface finishes after cutting
- | Optimal cutting shortens the further preparation process



PRECISION CLAMPING DEVICES


The Precision Series from QATM offers versatile clamping for both regular and irregularly shaped samples, designed for the most demanding and precise materialographic applications.

- | Handles Irregular and Complex Geometries
- | Varied and Versatile Product Range
- | Works with Rotating Device Integration (Qcut 150 A)
- | Horizontal and Vertical Adaptation Available



PRECISIE DOORSLIJPMACHINE QCUT 150 A

TECHNISCHE GEGEVENS

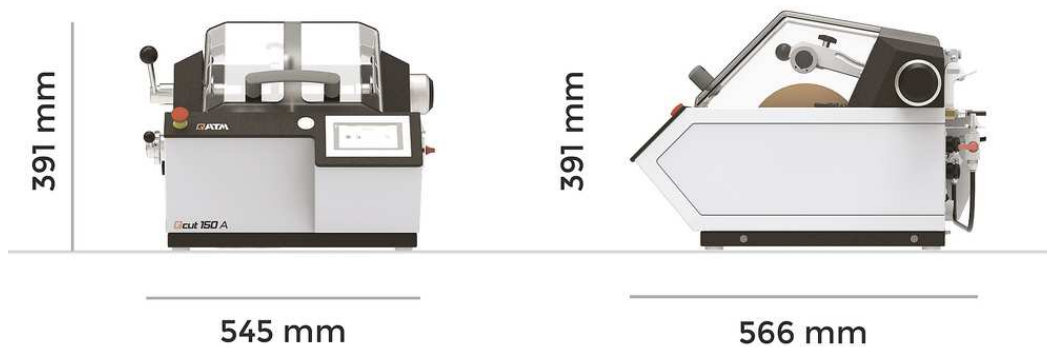

Slijpschijf	ø 75 - 203 mm/8"
Variabel toerental	150 - 3000 tpm
Max. monstergrootte	ø 40 mm
Schijfruimte	ø 12.7 mm
Variabel aanpasbare snij-kracht	0 - 10 N
Aanvoer	0,01 - 20 mm/s (automatisch)
Y-as (verticale beweging)	120 mm (manueel /automatisch)
Precision Y-axis	0.1 mm
Z-as manueel (parallel snede)	25 mm (manueel)
Precisie Z-as	0.005 mm
Motor vermogen	0.37 kW (S1)
Aansluitvermogen	0.75 kVA
B x H x D	545 x 391 x 566 mm
Gewicht	~ 45 kg (afhankelijk van uitvoering)

www.qatm.com/qcut-150-a

BESTELGEGEVENS

BASIS UITRUSTING

MACHINE (1 ITEM NOODZAKELIJK)

M1870020  Qcut 150 A (Brillant 210 A) Base machine 200-240V, 50/60Hz

M1870021  Qcut 150 A (Brillant 210 A) Base machine 100-120V, 50/60Hz

RECIRCULATION SYSTEM (1 ITEM REQUIRED)

A1870022 Qcut 150 A cutting wheel immersion cooling

A1870023 Qcut 150 A recirculation cooling system 15 Ltr.

A1870024 Qcut 150 A recirculation cooling system 45 Ltr., mobile

A1870025 Qcut 150 A recirculation cooling system 45 Ltr., for system lab

A1870026 Fresh water cooling system

OPTIONS (TO BE ORDERED WITH MACHINE)

A1870027 Dressing device

A1870028 Qcut 150 A preparation rotation device

OPTIES (KUNNEN AFZONDERLIJK BESTELD WORDEN)

Z1870031 Sample holder arm L standard

Z1870033 Sample holder arm L incl. cutting-force control

Z1870059 Rotation device type RD1

Z1870020 Table attachment

Z1870021 Qcut 150 A splash guard large angled

Z1870023 Parallel cutting guide

Z1870024 Angle cutting guide

Z1870022 Cutting guide fixed 90°

Z1870050 Flange set Ø 59 mm

Z1870051 Flange set Ø 99 mm

Z1870052	Flange set Ø 139 mm
Z1870035	Sample holder L standard Ø 40 mm open
Z1870037	Sample holder L standard Ø 40 mm closed
Z1870039	Sample holder L Flex Ø 40 mm
Z1870043	Sample holder L parallel clamping vice open
Z1870041	Sample holder L parallel clamping vice closed
Z1870045	Drop shape sample holder L Ø 40 mm
Z1870047	Universal clamping vice L
Z1870049	Matrix clamping vice L
Z1870054	Mini vice L
Z1870056	Double parallel vice L
Z1870058	Adaptor with ball joint L
Z1870064	Reduction ring Ø 25 mm
Z1870063	Reduction ring Ø 30 mm
Z1870071	Reduction ring Ø 40 mm
Z1870067	Reduction ring 1"
Z1870066	Reduction ring 1 1/4"
Z1870065	Reduction ring 1 1/2"
Z1870068	Manual chuck Ø 32 mm
Z1870069	Manual chuck Ø 44 mm

VACUUM CLAMPING VICES

Z1870028	Vacuum sample holder
01870288	Vacuum holder (sacrificial plate)
Z1304600	Vacuum pump 230V/50Hz for vacuum sample holder
Z1304634	Vacuum pump 110V/60Hz for vacuum sample holder

UNIVERSAL VICES AND ACCESSORIES

92005796	Longitudinal v-shoe
92005797	Traverse v-shoe
92005799	Longitudinal bridge

92005800	Traverse bridge
92005798	Standard clamping shoe
Z1350007	Height adapter 60mm Piccolo
Z1350003	Height adapter 74mm Mono
95003172	Clamping arm extension Piccolo 54 and Mono 61
Z1304591	Fixing screw M6 x 25
95003169	Fixing screw Piccolo M8 x 25
95003170	Fixing screw Piccolo M10 x 25
95003168	Fixing screw Mono M10 x 32
91001157	T-slot nut - stainless M6 x 8
91001161	T-slot nut - stainless M8 x 10
91001159	T-slot nut - stainless M8 x 12
91001158	T-slot nut - stainless M10 x 12