



MICRO / MACRO SCRATCH TESTER

QSCRATCH 20

The Qscratch 20 is a fully automatic macro/micro scratch tester for quality control of coatings, functional surfaces, and critical loads. The system combines controlled force application, precise electronic load control, integrated optical evaluation, and guided operation in QpixControl2. This enables reproducible scratch testing with progressive or constant load, including structured evaluation, documentation, and reporting.

With motorized X/Y/Z axes, single and multiple measurements, and an 8-position automatic tool changer, the Qscratch 20 is designed for laboratories that need reliable surface-related decisions every day. It supports scratch tests according to standards such as DIN EN ISO 20502, ASTM C1624, ISO 27307, and ASTM G171. Optional Vickers hardness testing and Qpix Inspect analysis modules make the system a flexible platform for coating evaluation, surface testing, and quality assurance.



MICRO / MACRO SCRATCH TESTER QSCRATCH 20

ПРЕИМУЩЕСТВА ОБОРУДОВАНИЯ

- | Fully automatic macro/micro scratch tester for coatings and surfaces
- | Single and multiple measurements for efficient test series
- | Electronic force control with a test force range of 0.5–200 N
- | Constant and progressive load for standards compliant scratch testing
- | 8-position automatic tool changer for test modules and lenses
- | Motorized X/Y/Z axes with precise positioning
- | QpixControl2 software with guided scratch workflow
- | Optional hardness upgrade and Qpix Inspect analysis functions

NOMINATED FOR THE BEST OF INDUSTRY AWARD

The Qscratch 20 has been nominated for the Best of Industry Award in the category "Best of Manufacturing & Design – Quality & Precision." Please support us with your vote!



MICRO / MACRO SCRATCH TESTER QSCRATCH 20

EXPERIENCE THE 3D MODEL IN THE REAL WORLD!

SHARE



CHOOSE PRODUCT:

AP-Model - Visit the page with your smartphone or scan the QR code under "View in Room" and experience the 3D model in the real world!

Qscratch 20 

MICRO / MACRO SCRATCH TESTER QSCRATCH 20

МЕТОДЫ ИЗМЕРЕНИЙ & ИСПЫТАТЕЛЬНАЯ НАГРУЗКА




METHOD: PROGRESSIVE LOAD

The initial and final loads, as well as the load rate, are defined. The load increases along the scratch length; visible events can be assigned as Lc points at their respective positions.

ISO: PFST - progressive force scratch test

ASTM: PL - progressive load

Lc: Load critical

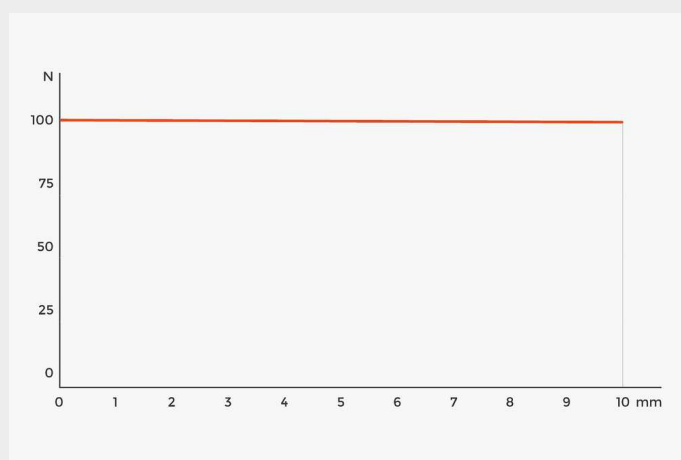

METHOD: CONSTANT LOAD

The normal force remains constant over the scratch length. This method is suitable for individual measurements, series comparisons, and recurring test conditions in a quality control setting.

ISO: CFST - constant force scratch test

ASTM: CL - constant load

Lc: Load critical

TEST METHODS

PROGRESSIVE OR CONSTANT
LOAD -
TAILORED TO THE TEST
REQUIREMENTS



STANDARDS & PROCEDURES

DIN EN ISO 20502 / ASTM C1624 -

Determination of adhesion and mechanical failure modes of hard ceramic coatings using a scratch test.

ISO 27307 - Evaluation of the adhesion and cohesion of thermally sprayed ceramic layers on a cross-section using a transverse scratch test.

ASTM G171 -

Determination of the scratch hardness of material surfaces using a defined diamond stylus, a constant normal force, and evaluation of the remaining scratch width.

FOR LABORATORIES THAT MAKE CRITICAL DECISION EVERY DAY

Qscratch 20 is a automated macro/ micro scratch tester for coatings and functional surfaces. The system combines controlled force application, integrated optical evaluation, and guided software operation and a report generation. It was developed specifically for laboratories that make critical surface-related decisions on a daily basis.

PRACTICAL QUALITY ASSURANCE & INTEGRATED WORKFLOWS

Designed for quality assurance: The system combines standardized scratch tests with high repeatability, efficient documentation, and intuitive operation. QpixControl2, the integrated workflow, flexible export options, and optional inspection functions turn it into a powerful, all-in-one solution rather than a simple stand-alone tester.

LONG-TERM INVESTMENT SECURITY THROUGH QATM EXPERTISE

QATM brings many years of experience in hardness testing, metallography, and global service to its scratch testing portfolio, thereby safeguarding your investment in the long term.

TEST PRINCIPLE DIN EN ISO 20502 | ASTM C1624

A DEFINED SCRATCH. A CONTROLLED LOAD. A CLEAR

EVALUATION.

The principle

In the scratch test, a diamond stylus is moved across the coated surface under a defined normal force. Visible events such as cracking, spalling, or penetration of the coating are localized along the scratch mark and assigned to the corresponding critical load. Thus, the scratch mark not only indicates that a failure is occurring; it also shows where it is occurring—and makes the load at the critical point traceable.



DEFINING AND ASSIGNING FAILURE EVENTS

Lc1 - Initial cracking / plastic deformation

Lc2 - Spalling

Lc3 - Penetration of the coating down to the substrate

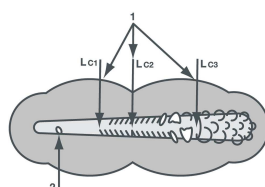


Figure: Failure Events

FROM SCRATCH TRACK TO TRACEABLE RESULT

CERAMIC COATINGS. EVALUATION OF ADHESION AND FAILURE MODES.

The following interactive section shows how typical failure events along the scratch track are analyzed. Click on the LC points to see how cracking, spalling, and coating penetration are associated with specific positions and loads.

SELECT A MATERIAL AND EXPLORE THE LC POINTS



**THIN, HARD CERAMIC
COATINGS (TYPICALLY 0.1 - 30
MM)**

- | Carbide, TiC, CrC
- | Nitride TiN, CrN, TiAlN
- | Oxide Al_2O_3
- | Diamond Diamond coating
- | Diamond-like carbon DLC

**THIN, HARD CERAMIC AND
TECHNICAL COATINGS
(TYPICALLY ≤ 20 MM)**

- | PVD TiN
- | PACVD DLC
- | Arc-discharge DLC
- | PVD/PACVD Cr-C
- | PVD Al_2O_3
- | Hard chrome

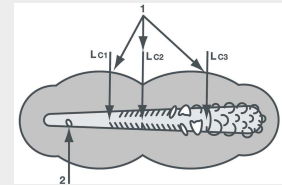



Figure: Failure
Events

MICRO / MACRO SCRATCH
TESTER QSCRATCH 20

РЕВОЛЮЦИОННАЯ ОПТИЧЕСКАЯ СИСТЕМА

The QATM-developed, in-house manufactured lens system sets new standards. As well as providing crystal clear image quality for hardness testing, Koehler illumination uses white LED light and motor-operated aperture shuttering to produce ideal contrast, even for high magnification images.

Experienced metallurgists agree the image quality provided by the Qscratch is comparable in all aspects with that of established sophisticated microscopes. The up-to-date concept and new lenses in the optic system enable the device to completely meet even the strictest physical 'test system definition' requirements in compliance with DIN EN ISO6507-1/2:2018.

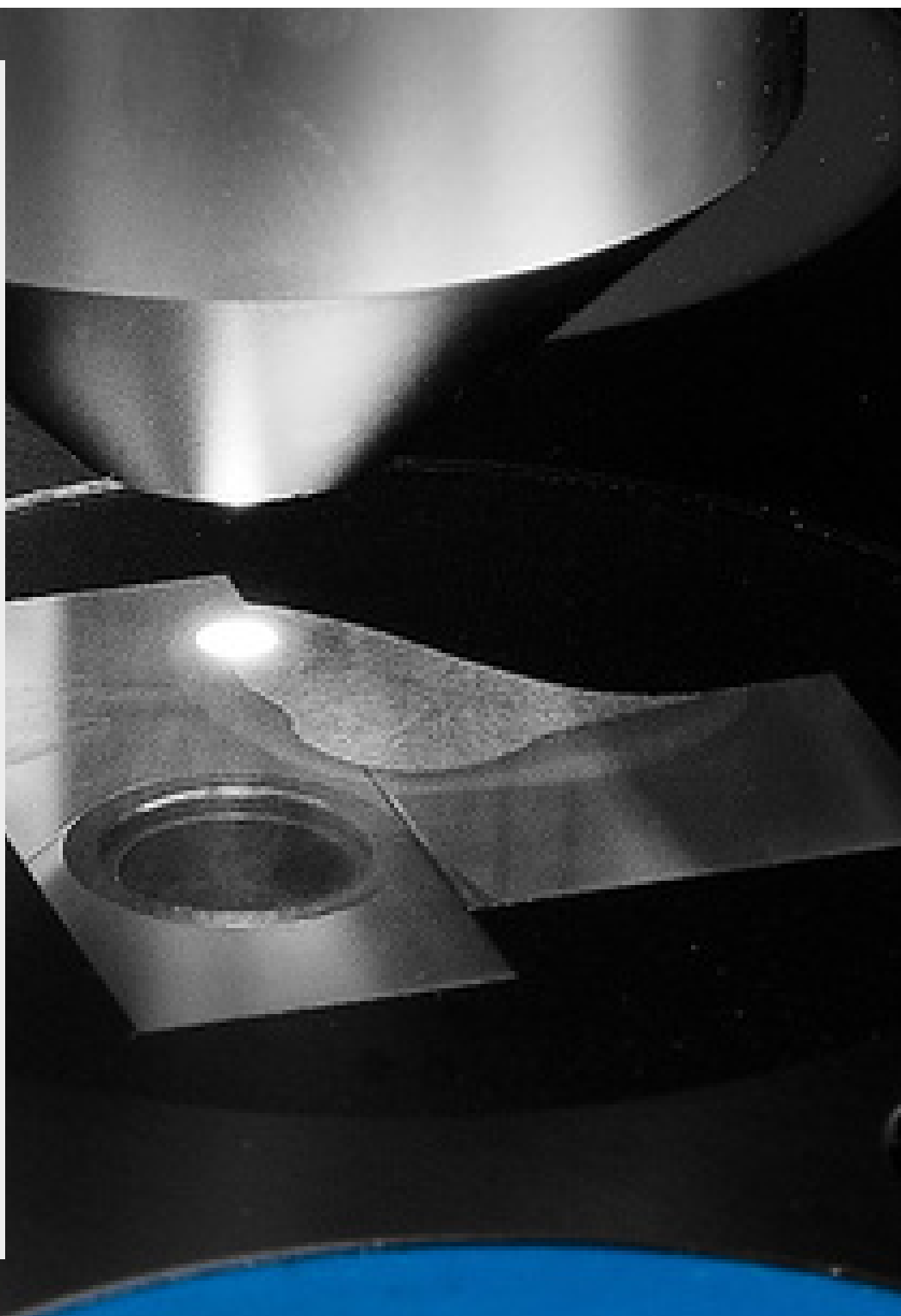

MICRO / MACRO SCRATCH TESTER QSCRATCH 20

INNOVATIVE OPERATION



ПОДСВЕТКА ТЕСТОВОГО ПРОСТРАНСТВА

Все приборы оснащены новым светодиодным освещением рабочего пространства: упрощенное позиционирование образцов для единичных измерений.

ИНДИКАЦИЯ СОСТОЯНИЯ
С ПОДСВЕТКОЙ

ОСВЕЩАЕТ ПРОСТРАНСТВО


The illuminated QATM logo displays the current device status at a glance. The range of flash intervals indicates whether the device is operating automatically or is free to be used for new tasks for staff all around the lab. In addition, the LED test space lighting, installed as standard, ensures that samples and sample holders are positioned correctly.

On the Qscratch 20, it also ensures uniform light intensity for sample imaging.



MICRO / MACRO SCRATCH TESTER QSCRATCH 20

ТЕХНОЛОГИИ ПИОНЕРА - УНИКАЛЬНАЯ РЕАЛИЗАЦИЯ

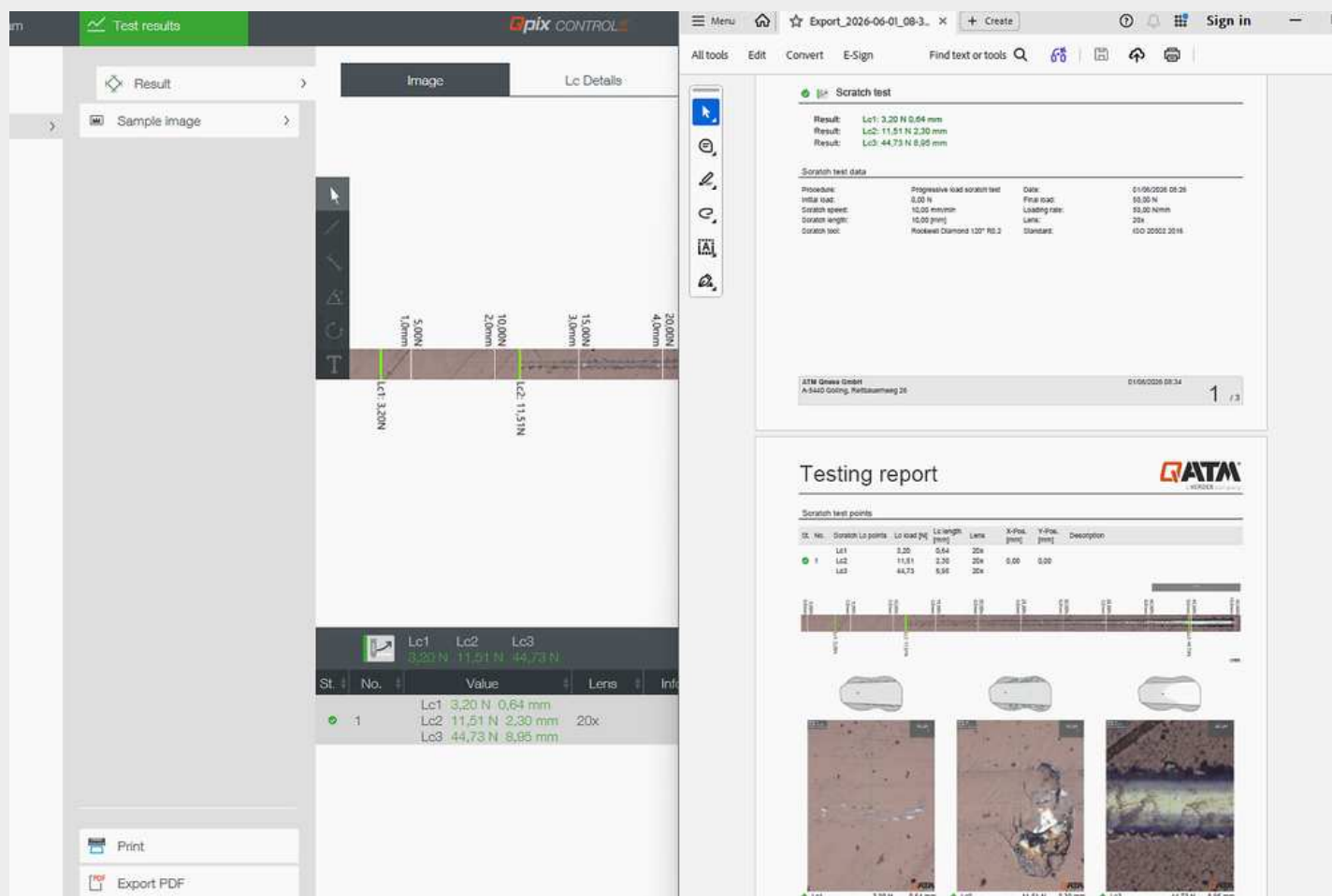


EVALUABLE AND COMPARABLE

The selected loading logic is linked to the scratch mark, position, and result. Methods can be applied to sample series, coating variants, and process comparisons.




CHECKING THE DIAMOND STYLUS

The diamond tip can be visually inspected directly within the system. This allows the condition of the diamond stylus to be checked before or after the test - an important step for ensuring reproducible results and traceable test conditions.

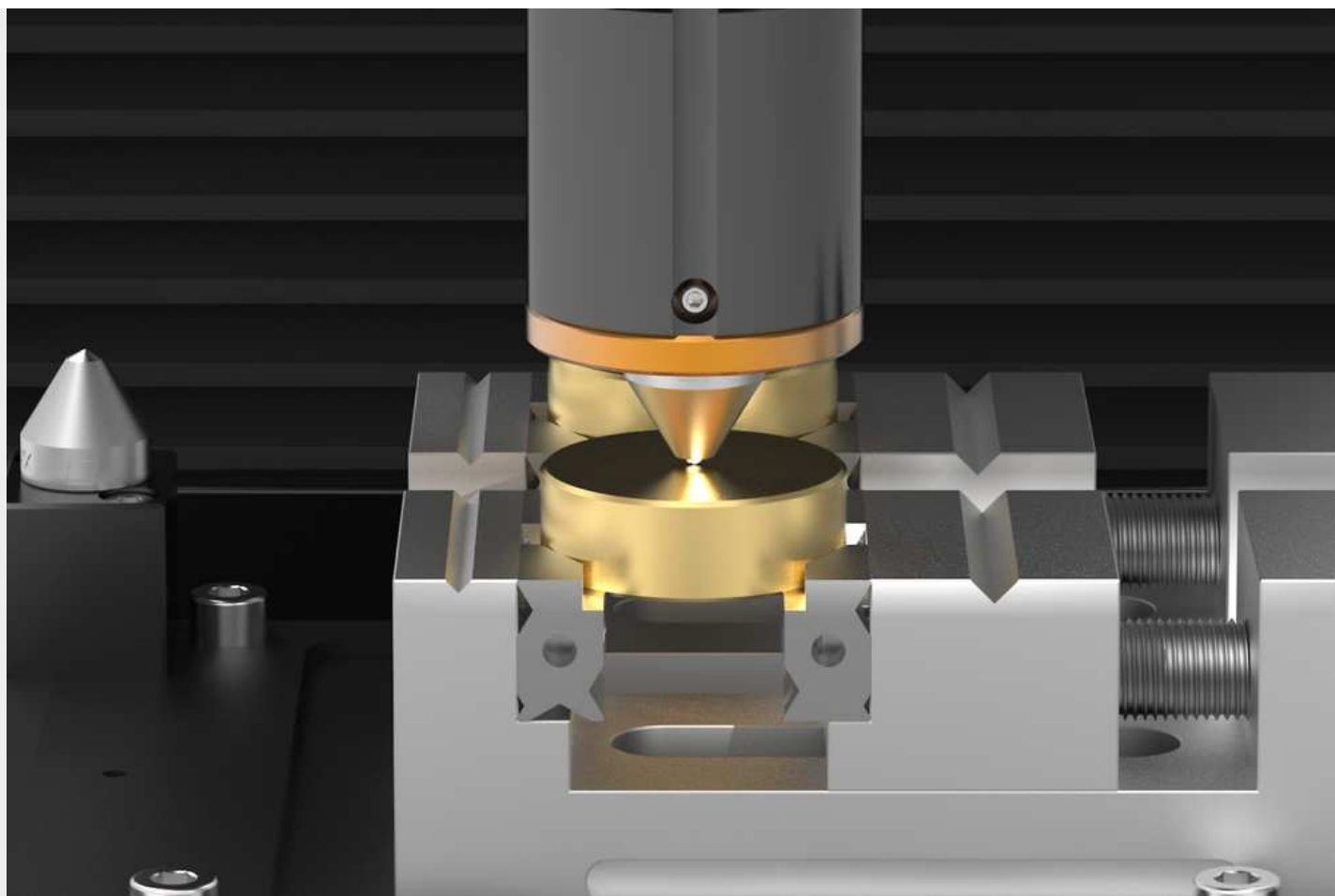


UNIVERSAL SAMPLE HOLDER

Adjustable clamping jaws allow for the secure holding of various specimen geometries - from flat and angular test pieces to round or very small specimens.

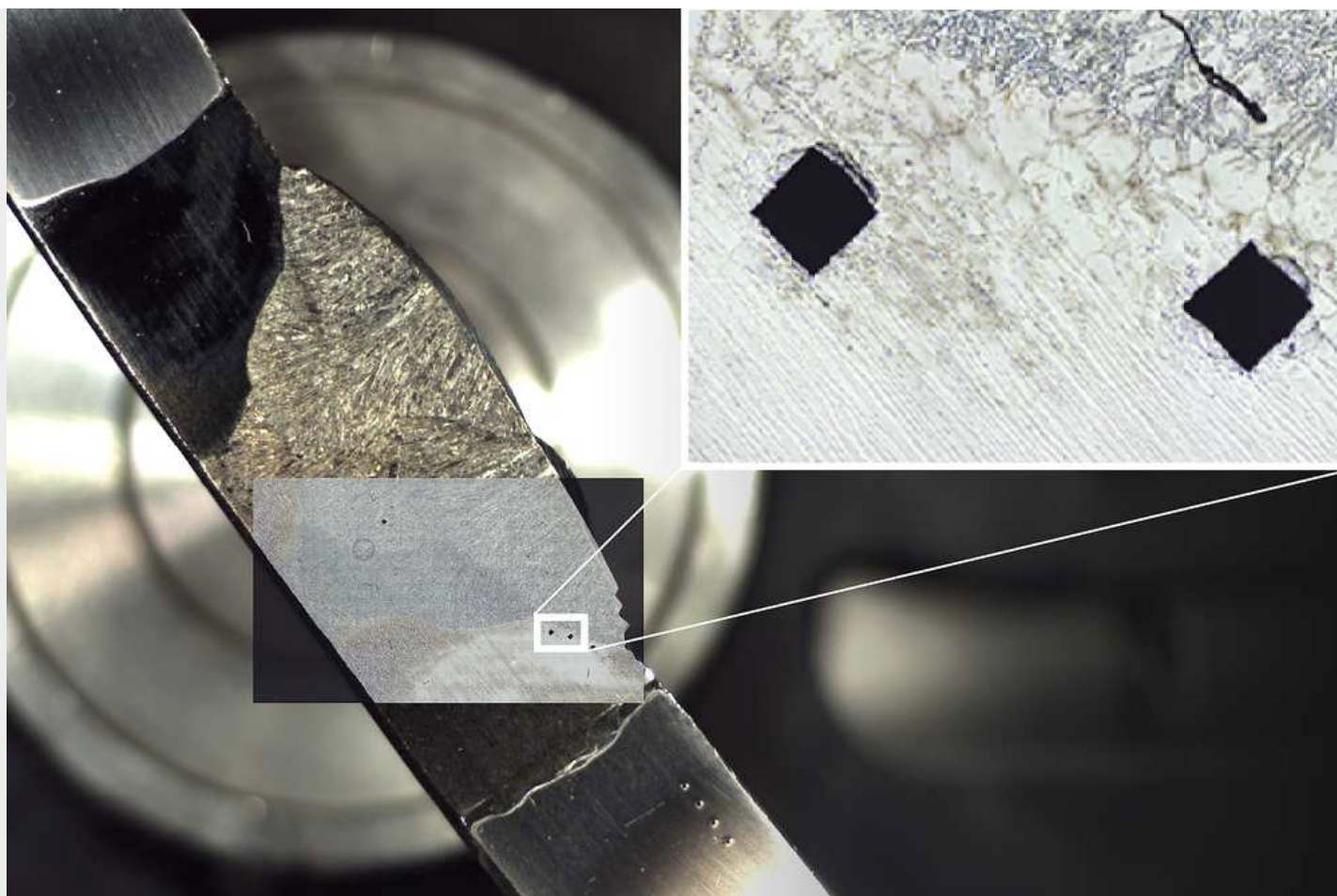


HIGH-RESOLUTION SAMPLE IMAGE (HRI)

If high-quality images of larger areas are required (e.g. for weld seam measurements), the area can be scanned using the HRI function. The Qpix Control 2 software automatically combines the individual images into one large overall image.

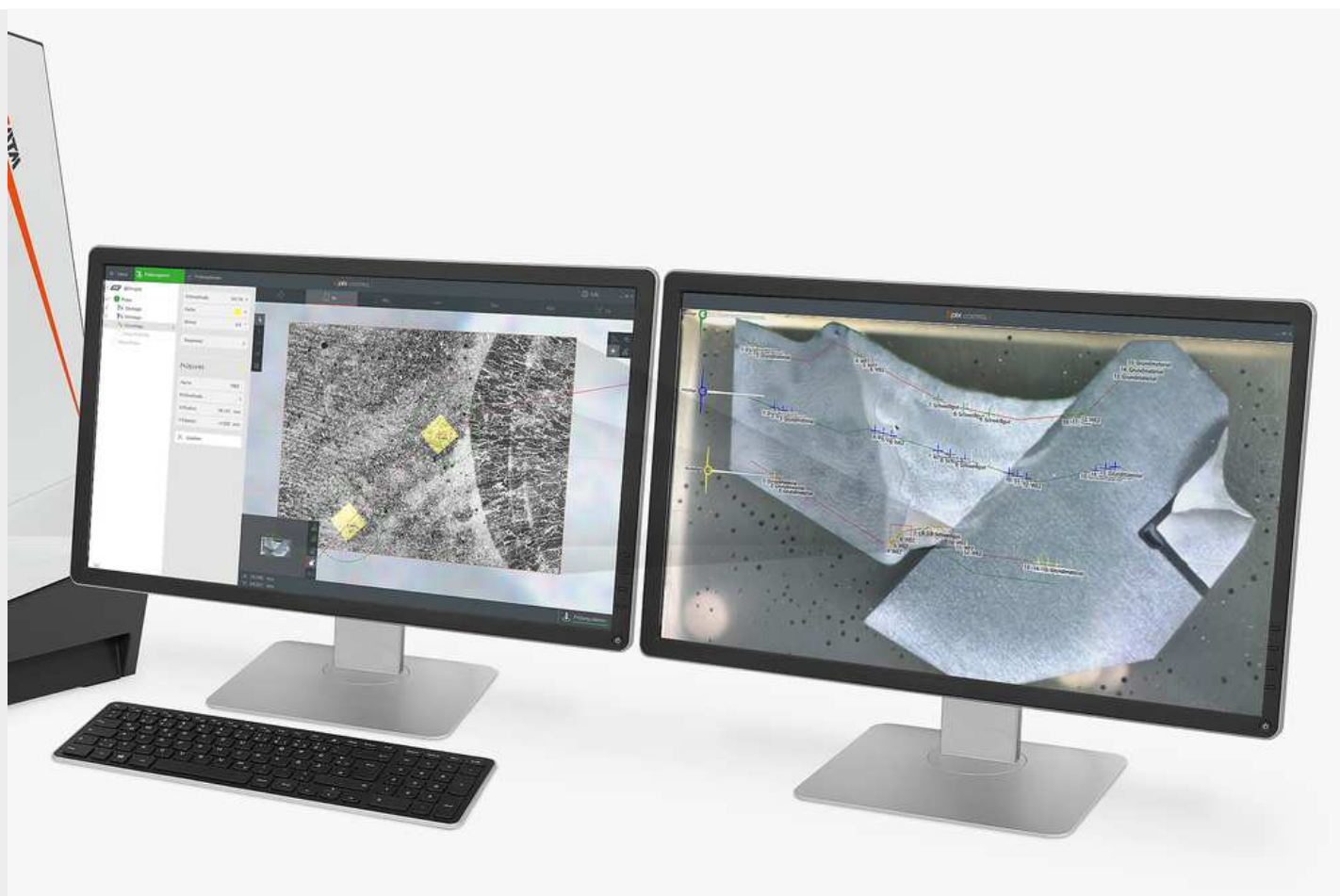

ОБЗОРНАЯ КАМЕРА НАБЛЮДЕНИЯ ОБРАЗЦА

It's no coincidence that most QATM customers choose the 'A+' version with a built-in sample image camera. In a few seconds the image of the sample is shot with the additional camera (field of view 52 x 39 mm). The image provides excellent navigational support within the software, particularly in combination with DOUBLE-VIEW TECHNOLOGY, and aids enhanced documentation in the automatically compiled test report.

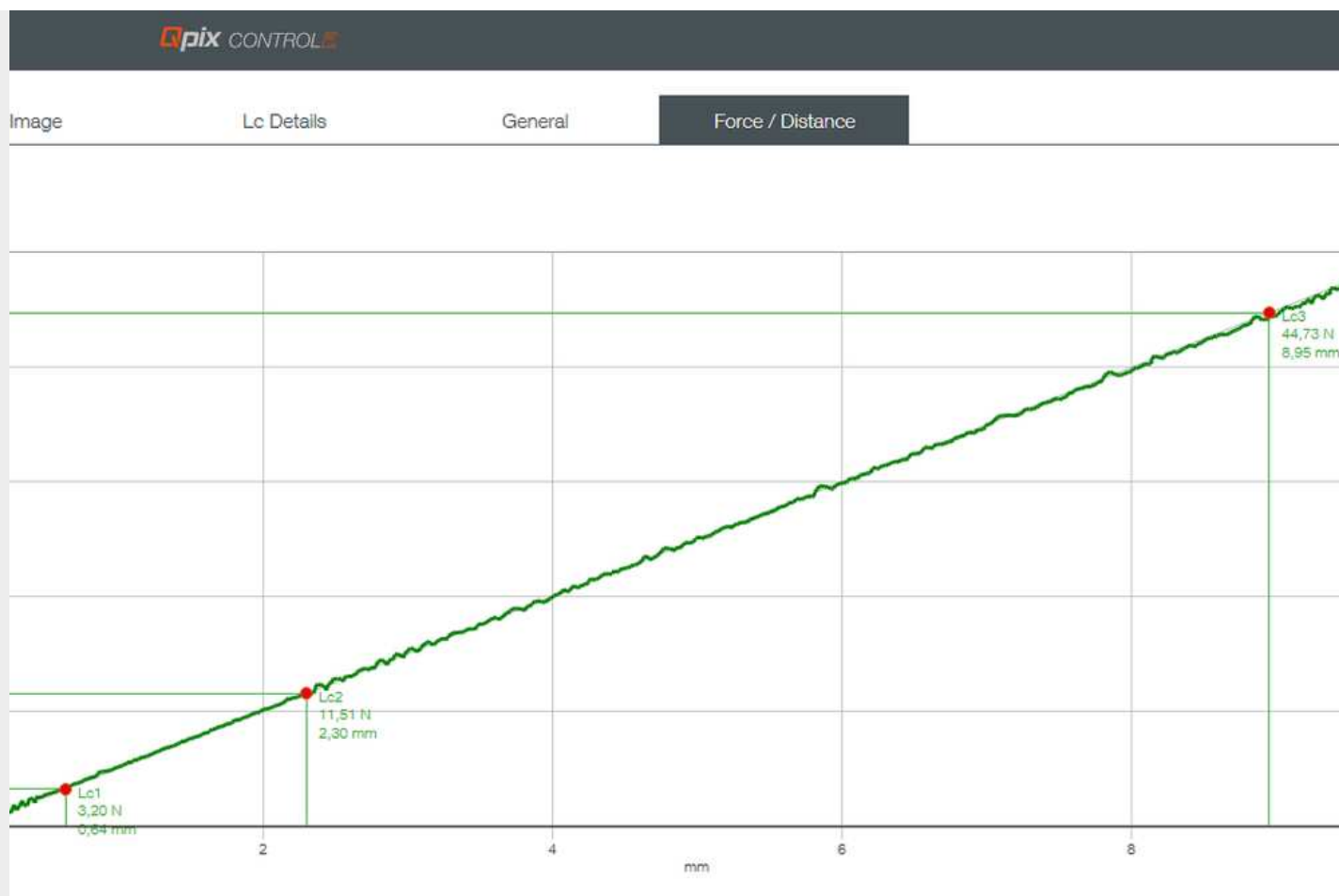

PRECISE DISTANCE / FORCE APPLICATION


Precise electronic load control ensures that the test load—as a critical test parameter—remains controlled, reproducible, and traceable at all times.

Qpix CONTROL²

SCRATCH WORKFLOW

EASY TESTING, EFFICIENT ANALYSIS

QpixControl2 is more than just a user interface: The software combines a patented 3D user concept with a seamless scratch workflow. Tests can be easily created, scratch marks visually positioned, single or multiple measurements processed automatically, and relevant Lc events evaluated in a structured manner. From test definition through image capture and evaluation to data management, export, and standardized reporting, all steps remain integrated within an intuitive software environment. This reduces operational effort, supports reproducible processes, and enables the efficient use of scratch testing in production-oriented quality assurance.

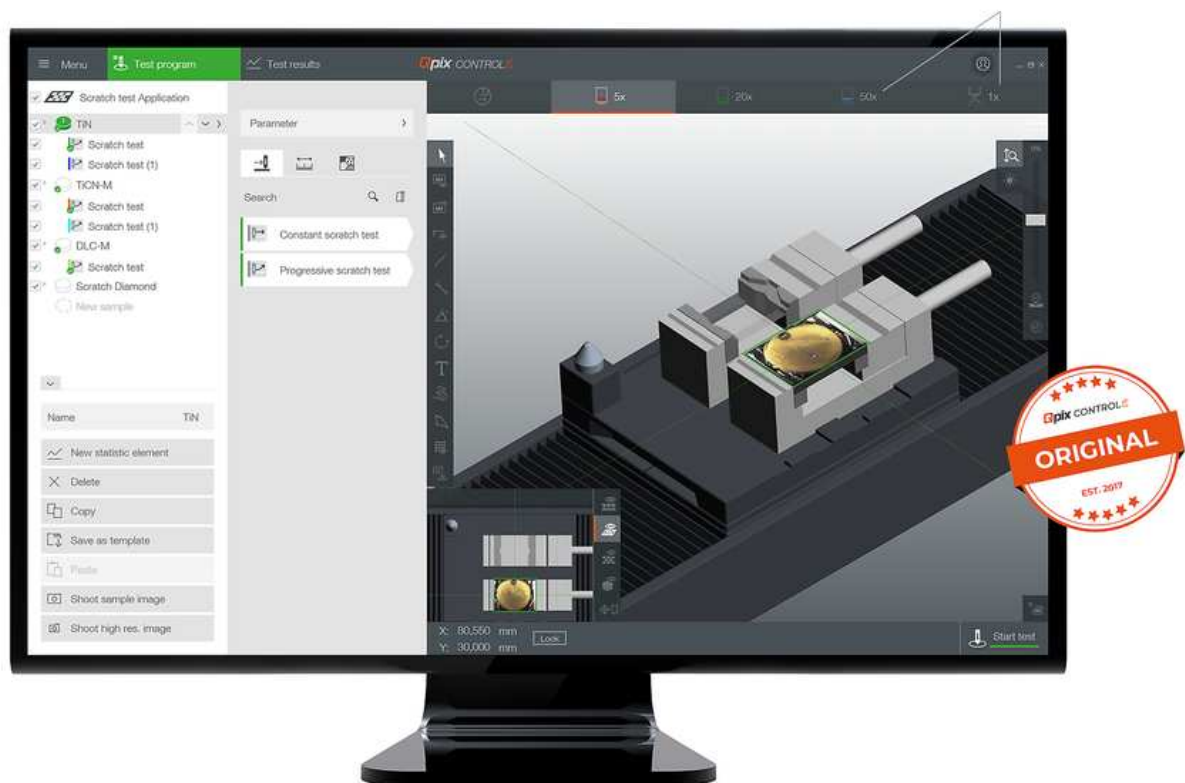


FROM SETUP TO REPORT 3 ШАГА ДО РЕЗУЛЬТАТА




1. LOAD SAMPLE

The sample holder is automatically moved to the correct height, and the sample image is automatically captured.

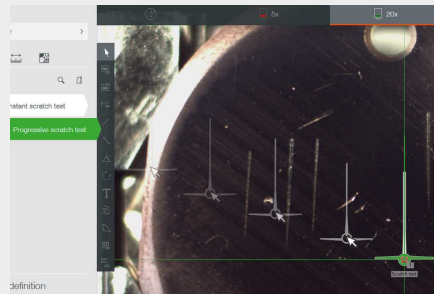

2. POSITIONING


Drag and drop the starting point of the scratch test directly to the desired position.

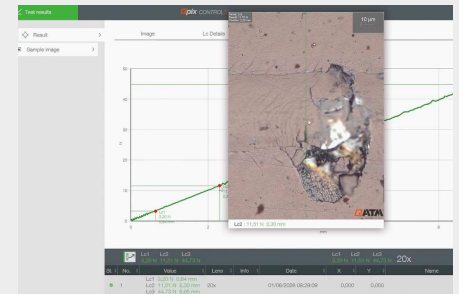

3. ЗАПУСК ЦИКЛА ИЗМЕРЕНИЙ

The test sequence is fully automatically and in accordance with the Scratch test standards. Subsequently, the failure events can be evaluated and assigned to the respective Lc points.


ПРОГРАММНОЕ ОБЕСПЕЧЕНИЕ INCLUDED FEATURES

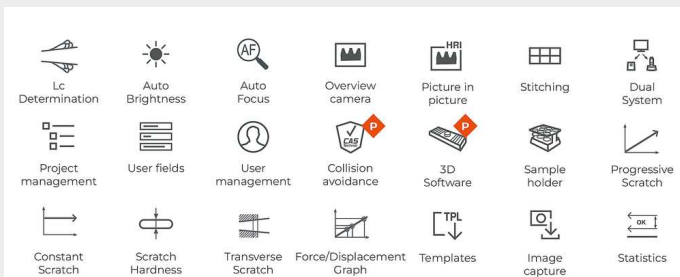


ПРОГРАММНОЕ ОБЕСПЕЧЕНИЕ ОПЦИИ



MICROSCOPY & ANALYSIS

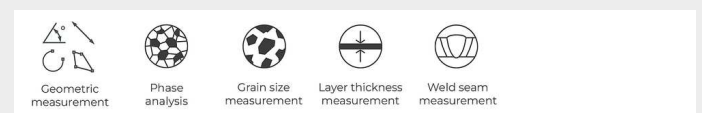


QPIX CONTROL2 INNOVATIVE SOFTWARE FUNCTIONS

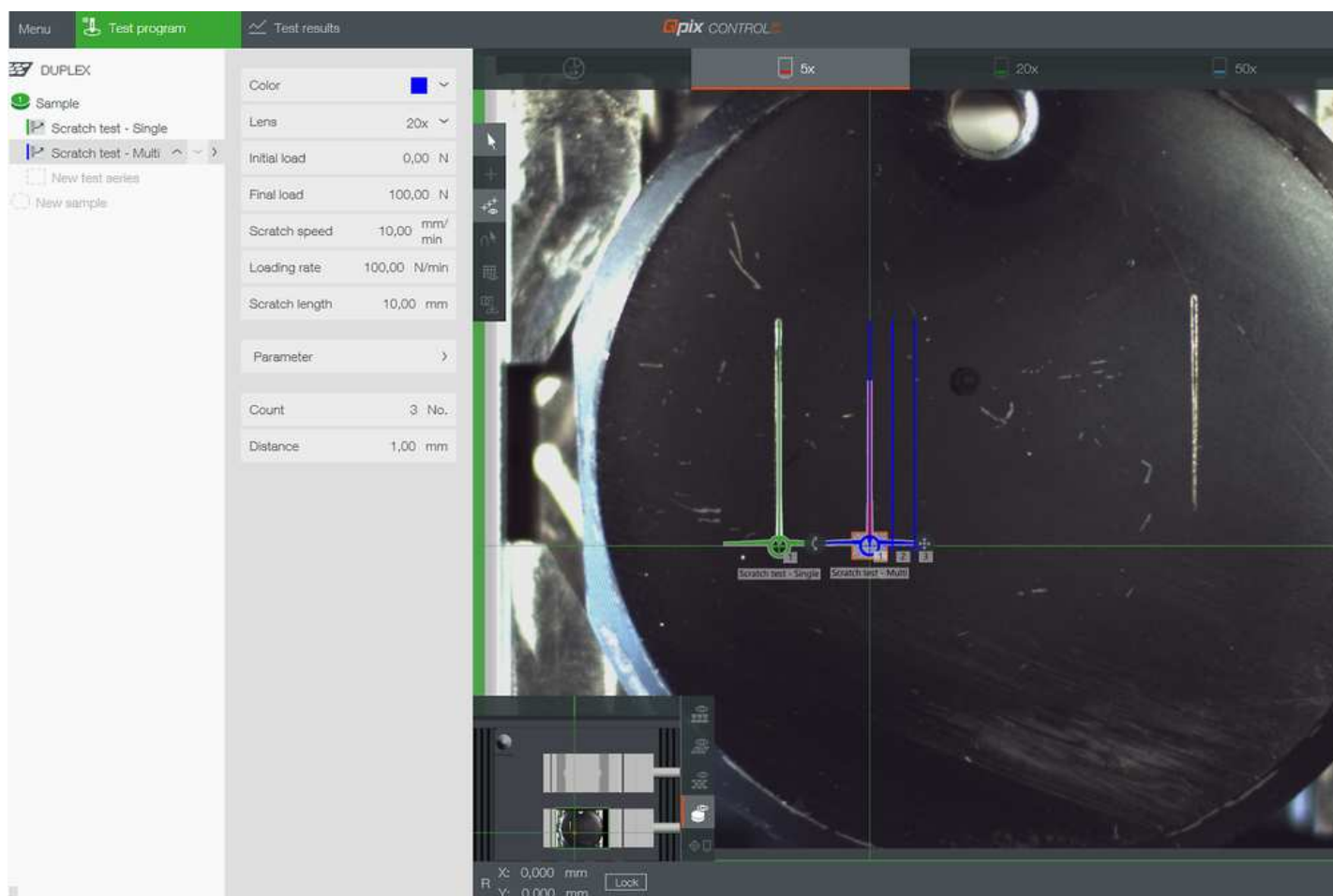


SINGLE AND MULTIPLE INSPECTIONS

Single or multiple measurements: Scratch marks can be positioned directly in the software, created with a defined offset, and processed automatically. The Lc results are evaluated for each test as well as as a series mean.

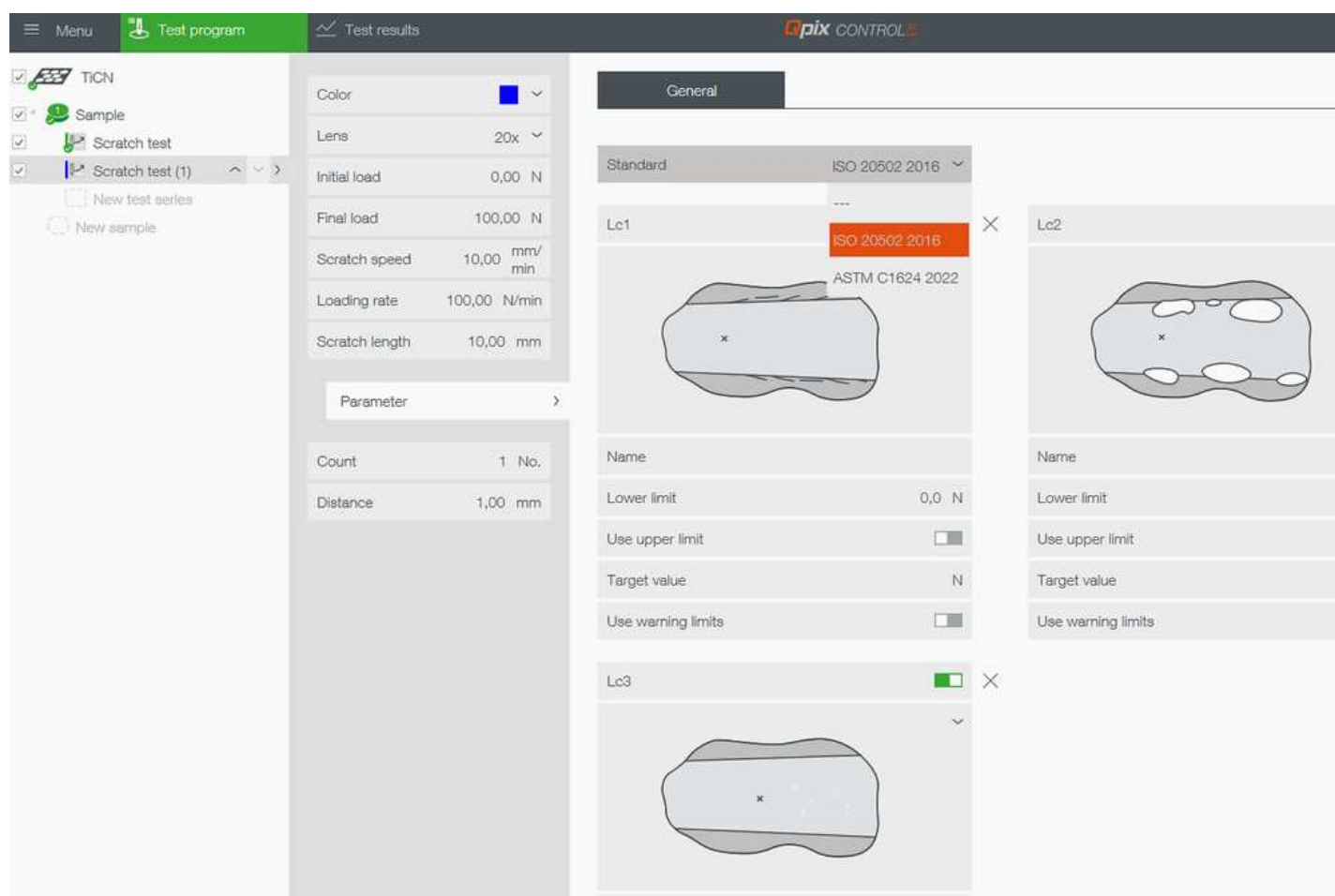


A WIDE RANGE OF SETTING OPTIONS

Load definition, crack length, load rate, speed, and lens selection are managed directly within the method. For evaluation, standardbased templates or custom templates can be used and supplemented with limit and target values.

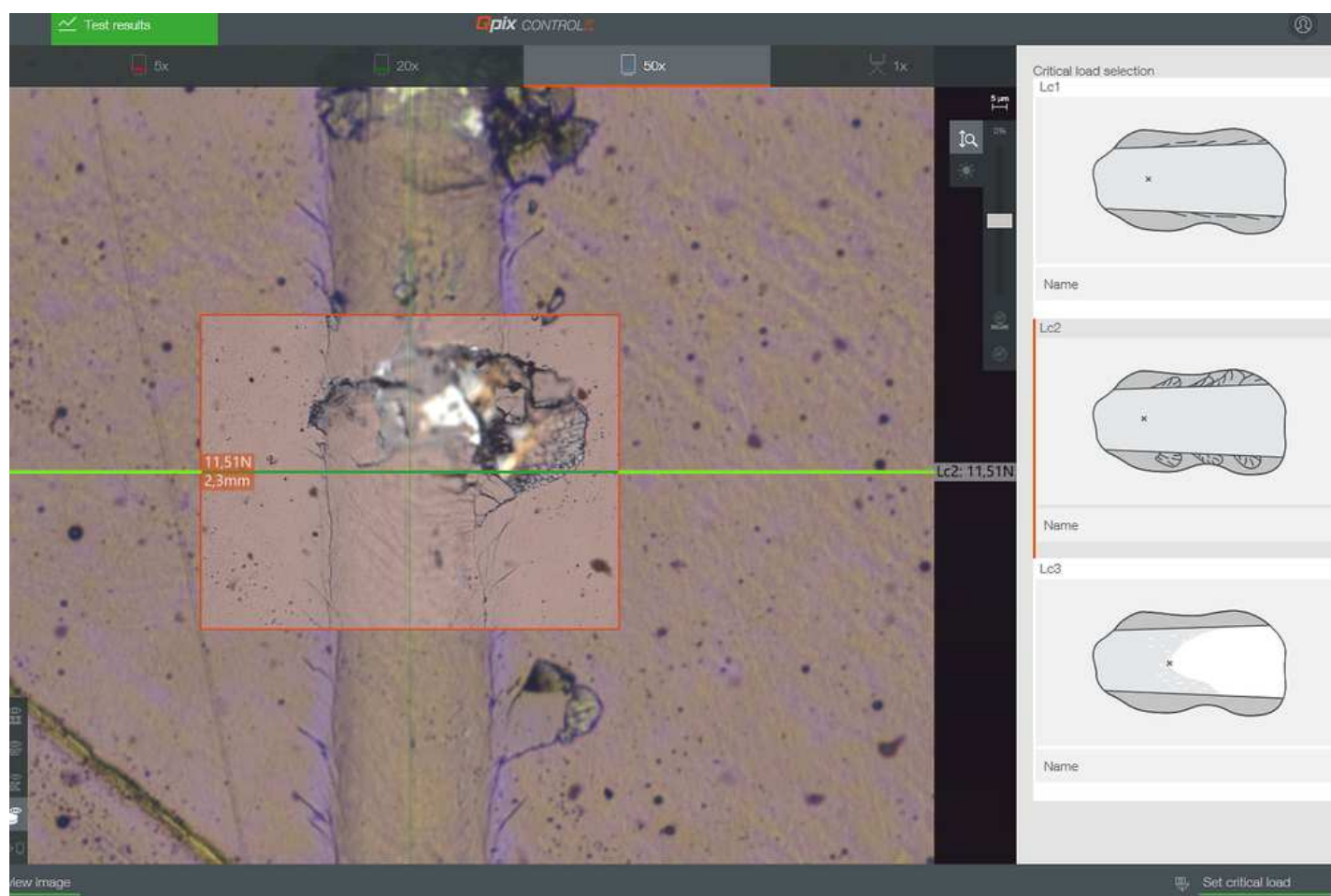

GUIDED LC DETERMINATION

Graphical templates support the consistent classification of typical failure events. This facilitates standard-based procedures and reduces room for interpretation.

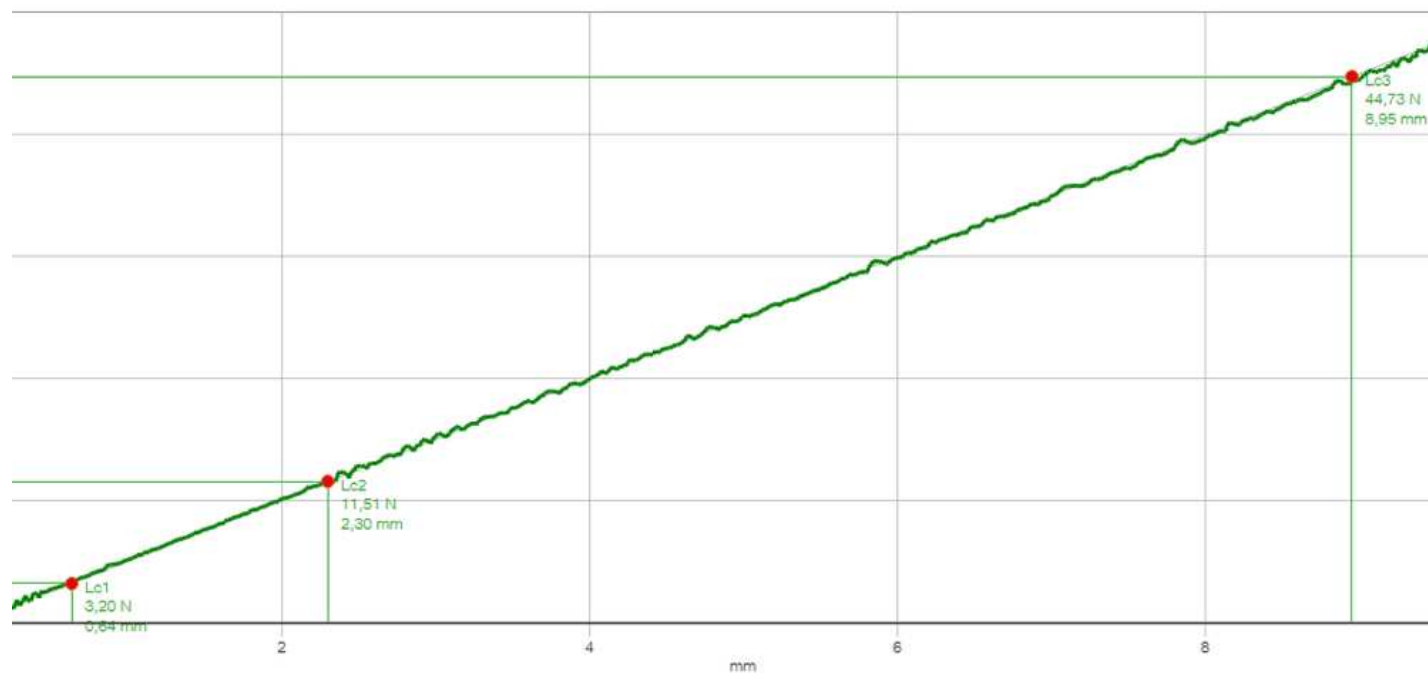


FORCE/DISTANCE GRAPH

The force-distance curve documents the applied load along the scratch track. Anomalies in the curve can be compared with visible failure events and used for Lc evaluation


results **pix CONTROL**

esult > Image **Lc Details** General Force / Distance


ple image >

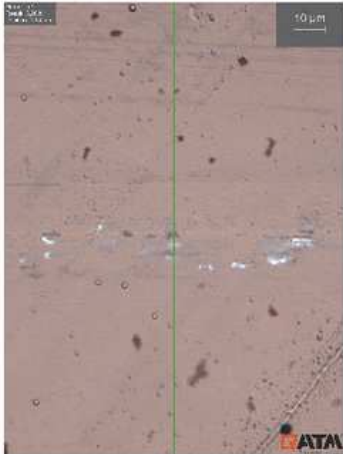

✓ Lc1 : 3,20 N 0,64 mm
Longitudinal cracks at the edges of the scratch track

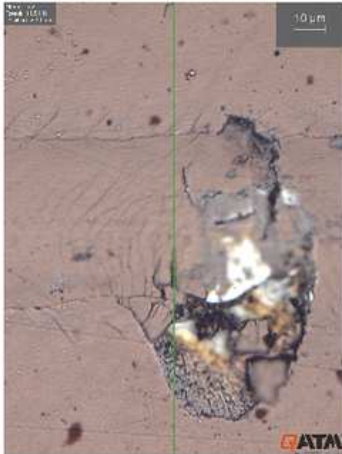

✓ Lc2 : 11,51 N 2,30 mm
Cohesive spallation along the edge of the scratch track

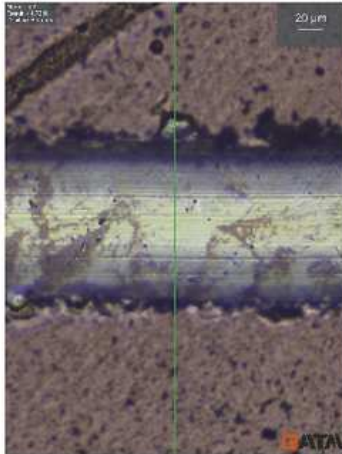

✓ Lc3 : 44,73 N 8,95 mm
Continuous ductile coating perforation



Lc1 Lc2 Lc3
3,20 N 11,51 N 44,73 N

Lc1 Lc2 Lc3
3,20 N 11,51 N 44,73 N 20x

LC DETAILS


Relevant failure events along the scratch track can be viewed in detail in the image and assigned to the corresponding Lc details. This transforms the visual assessment into a verifiable test result

Scratch test data

Procedure:	Progressive load scratch test	Date:	01/06/2026 08:26
Initial load:	0,00 N	Final load:	50,00 N
Scratch speed:	10,00 mm/min	Loading rate:	50,00 N/min
Scratch length:	10,00 [mm]	Lens:	20x
Scratch tool:	Rockwell Diamond 120° R0.2	Standard:	ISO 20502 2016

Scratch test points

St.	No.	Scratch Lc points	Lc load [N]	Lc length [mm]	Lens	X-Pos. [mm]	Y-Pos. [mm]	Description
✓	1	Lc1	3,20	0,64	20x	0,00	0,00	
		Lc2	11,51	2,30	20x			
		Lc3	44,73	8,95	20x			

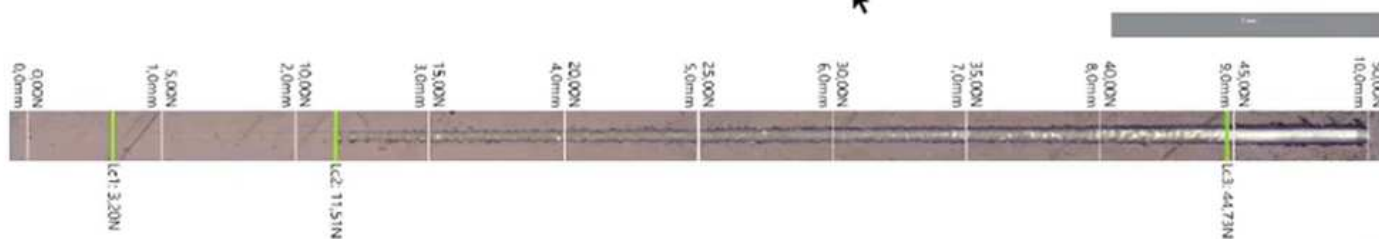

REPORT




Test parameters, scratch track, Lc values and evaluation results are compiled into a test report. This simplifies documentation, the comparison of test series, and the dissemination of results within the quality process.

ONE PLATFORM

MORE THAN JUST SCRATCH TESTING

**QSCRATCH 20 / QNESS 60 A+
EVO**

Qscratch 20 is part of a modular QATM platform. In addition to scratch testing, hardness testing methods, analysis functions, and other coating tests can be optionally integrated. This reduces the need to switch systems, consolidates results in a familiar software environment, and ensures long-term investment security.

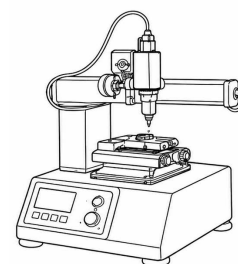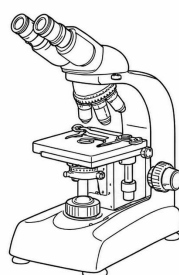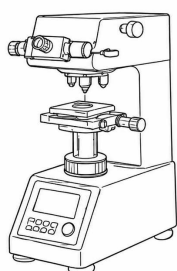

ИЗМЕРЕНИЕ ТВЕРДОСТИ

Optional hardness upgrade for Vickers, Knoop, Brinell, and Rockwell - integrated into the QpixControl2 software environment.

ANALYSIS

Optional analysis functions for coating thickness, phases, grain size, and weld seams - for more information from a single system.

SCRATCH TESTING

Controlled scratch testing with progressive or constant load - for evaluating coatings, surfaces, and critical loads.


ИЗМЕРЕНИЕ ТВЕРДОСТИ



Vickers

DIN EN ISO 6507, ASTM E-384, ASTM E92

HV 0.00025*	HV 0.0005*	HV 0.001	HV 0.002			
HV 0.003	HV 0.005	HV 0.01	HV 0.02	HV 0.025*		
HV 0.05	HV 0.1	HV 0.2	HV 0.3	HV 0.5	HV 1	HV 2
HV 2,5	HV 3	HV 5	HV 10	HV 20	HV 30	HV 50

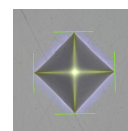


HV 60*



Knoop

DIN EN ISO 4545, ASTM E-384, ASTM E92


HK 0.001	HK 0.002	HK 0.005	HK 0.01	HK 0.015	
HK 0.02	HK 0.025	HK 0.05	HK 0.1	HK 0.2	HK 0.3
HK 0.5	HK 1	HK 2			

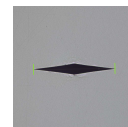

Brinell

DIN EN ISO 6506, ASTM E-10


HBW 1/1	HBW 1/2.5	HBW 1/5	HBW 1/10	HBW 1/30
HBW 2.5/6.5	HBW 2.5/31.25	HBW 2.5/62.5		
HBW 5/25	HBW 5/62.5			

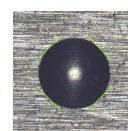

Rockwell

DIN EN ISO 6508, ASTM E-18

HRA	HRF	HR15-N/T	HR30-N/T	HR45-N/T
-----	-----	----------	----------	----------

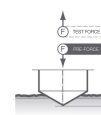

Перевод в другие шкалы согласно


DIN EN ISO 18265, DIN EN ISO 50150, ASTM E-140

* нестандартизированные измерения

COATING INSPECTION MADE EASY

QPIX INSPECT FOR CRATER-GRINDING AND ROCKWELL ADHESION TESTING

The optional Qpix Inspect modules expand the Qscratch 20 with practical functions for coating evaluation. Crater-grinding supports microscopic coating thickness measurement in accordance with DIN EN ISO 26423, while the Rockwell adhesion test enables qualitative adhesion evaluation of ceramic coatings in accordance with DIN EN ISO 26443. This allows additional information on coating quality to be captured, evaluated, and documented directly within the familiar QpixControl2 software environment.

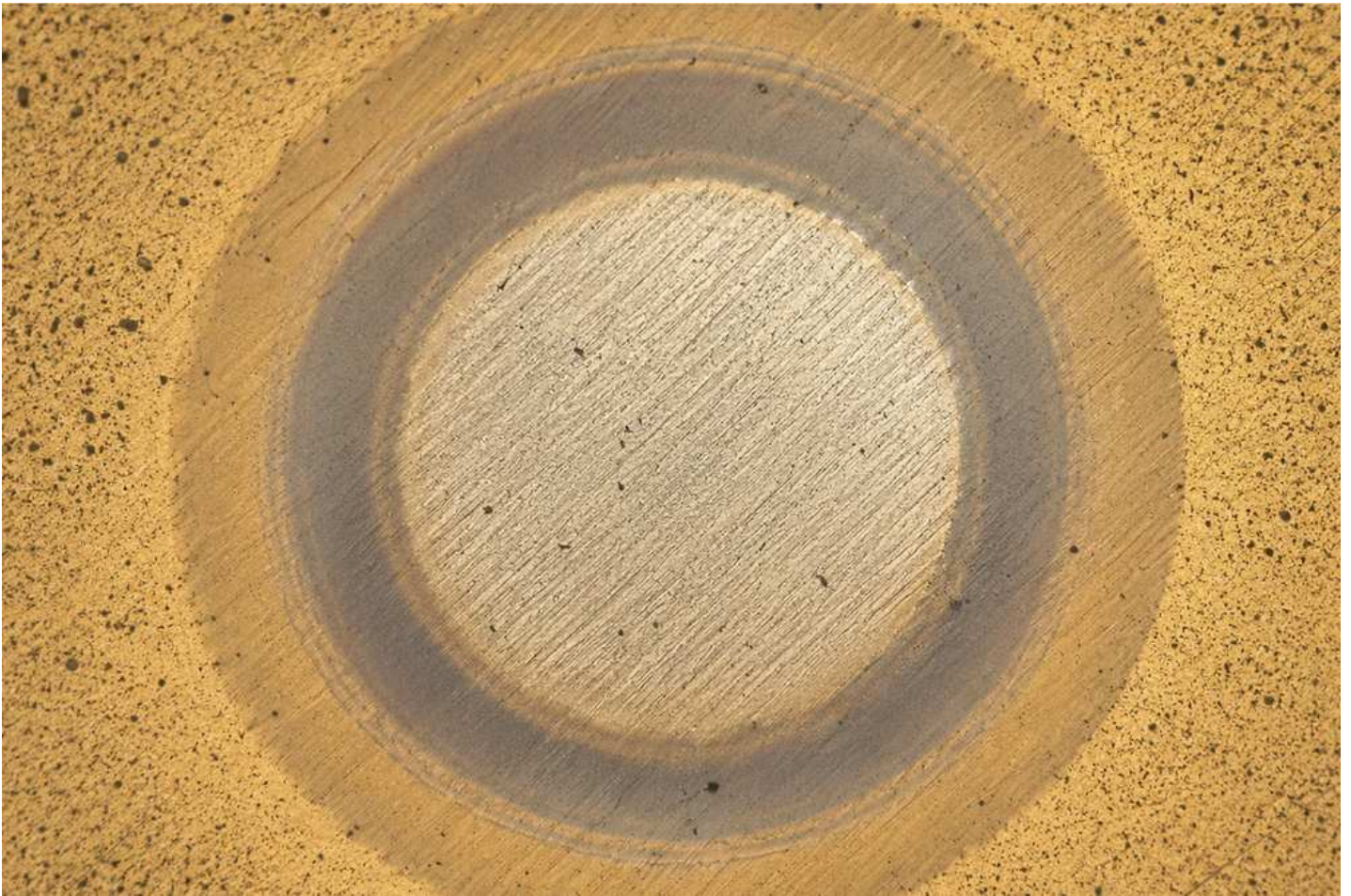

CRATER-GRINDING COATING THICKNESS MEASUREMENT

Crater-grinding method for microscopic coating thickness determination as supplementary information for coating evaluation in accordance with DIN EN ISO 26423.




ROCKWELL INDENTATION TEST

Rockwell indentation test for qualitative adhesion evaluation of ceramic coatings in accordance with DIN EN ISO 26443.

СТРУКТУРНЫЙ АНАЛИЗ СДЕЛАЛИ ДОСТУПНЫМ


МОДУЛИ ПРОГРАММНОГО ОБЕСПЕЧЕНИЯ QPIX INSPECT

The intuitive and user-friendly Qpix INSPECT software functionality provides a comprehensive toolbox for microscopic evaluations and result documentation. The multifunctional software can be customized for customer-specific measuring tasks and complemented with add-on modules.

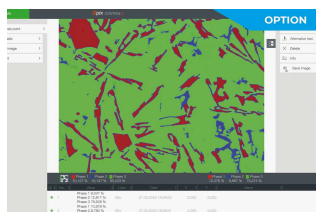


INSPECT ФАЗОВЫЙ АНАЛИЗ

- | Автоматическое определение размеров объектов изображения
- | Оценка фазовых фракций согласно ISO 9042 и ASTM E562
- | Предоставляет результаты аналитики в виде процентных долей поверхности или номинальных значений поверхности в виде таблиц или диаграмм

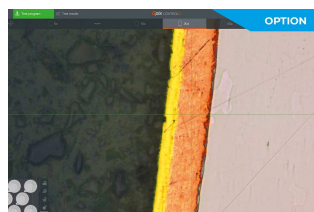


INSPECT ИЗМЕРЕНИЕ ТОЛЩИНЫ СЛОЕВ И ПОКРЫТИЙ

- | Измерение толщины слоев согласно DIN EN ISO 1463
- | Полуавтоматический замер горизонтальных, вертикальных и радиальных слоев
- | Предоставление толщины слоя в виде статистических линейных значений в виде таблиц или диаграмм

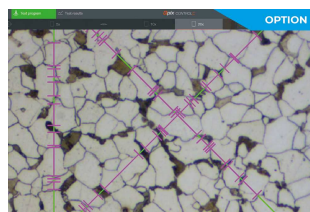

INSPECT ОПРЕДЕЛЕНИЕ РАЗМЕРА ЗЕРНА

- | Размер частиц определяется в соответствии с DIN EN ISO 643 и ASTM E112 методом линейного или кругового сечения.
- | Результаты анализа представляются в виде таблиц или диаграмм
- | Документирование статистических характеристик размера частиц и длин сегментов, пересекающих частицы

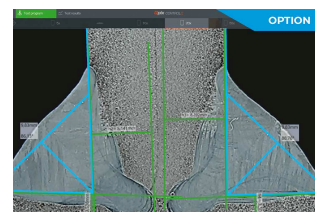


INSPECT WELD SEAM MEASUREMENT

- | Standardised measurement and evaluation of weld seams
- | Prefabricated templates with all relevant measuring tools such as throat thickness, weld reinforcement, penetration depth, etc.
- | Automatic good/bad evaluation and report generation.

IOT - INTERNET OF THINGS

THE PLATFORM FOR REMOTE ACCESS TO YOUR DEVICES

All QATM hardness testers with QpixControl2 and QpixT2 software seamlessly integrate into the Verder Scientific IoT platform, providing enhanced functionality and seamless connectivity.

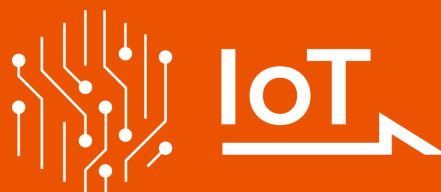



- | **Real-time Monitoring:** Monitor your machinery in real time, from anywhere in the world. This data-driven approach empowers you to make informed decisions with ease.
- | **Live Notifications:** Be ahead of the curve with immediate alerts and updates. Real-time notifications ensure you stay informed about your equipment's performance, leading to proactive maintenance.
- | **Effortless Backup:** Simplify your data protection. Whether you need to back up a single device or an entire fleet, our platform streamlines the process, minimizing downtime and data loss.
- | **Automatic & Free Software Updates:** Bid farewell to manual updates! Verder Scientific IoT ensures your customers' machines are consistently equipped with the latest software, optimizing performance and reliability.



ИНДУСТРИЯ 4.0



ДЛЯ ВОЗМОЖНОСТИ ИНТЕГРАЦИИ В СЕТЬ ЗАВТРА

Qconnect - это интерфейс в программном обеспечении QATM Qpix Control2, предоставляющий заказчикам полный набор возможностей обеспечения связи между приборами - от серийного производства, открытых XML-интерфейсов (двунаправленных) и предварительно определенных подключаемых решений, таких как QDAS Plug-In + до индивидуальных решений для заказчиков, полностью реализуемых компанией QATM. У нас есть профессиональное решение для каждого применения.

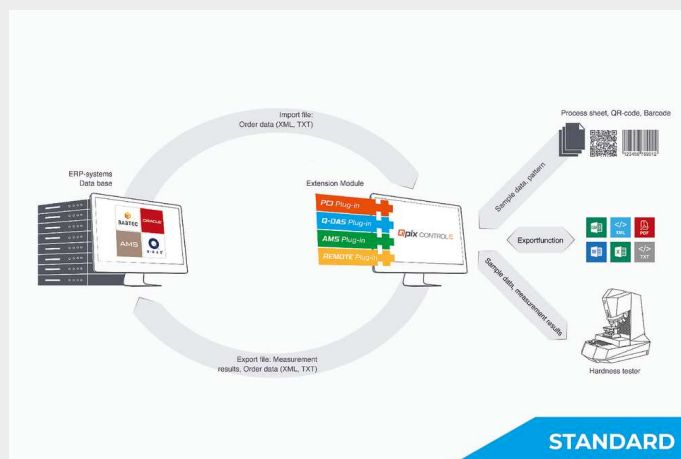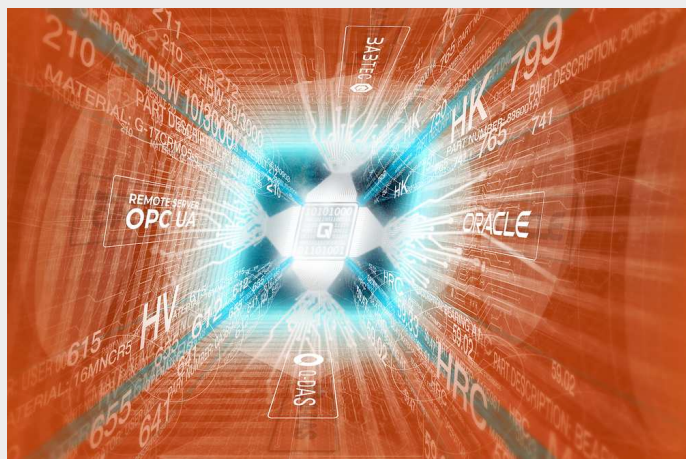

APPROVED QUALITY

DATA EXCHANGE & MONITORING

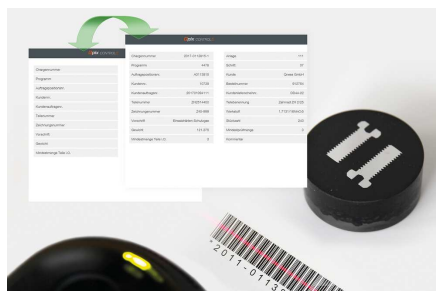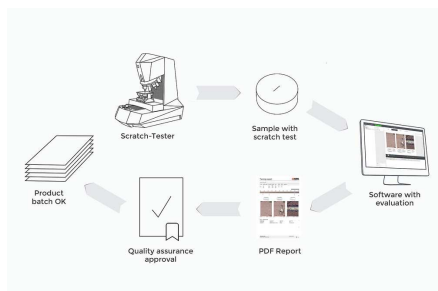

REPORT & AUTO-EXPORT


Test results, Lc results, scratch track information, and test parameters can be documented in a structured format and exported. Standardized test reports and flexible export options facilitate the transfer of results to quality assurance, documentation, or higher-level systems.

ШТРИХКОД/QR КОД/DMC СЧИТЫВАТЕЛЬ

Программные платформы Qpix поддерживают считыватели штрих-кода и QR-кода. Будь то простая вставка заголовочных файлов (последовательных), управление полной интеграцией считывающих устройств для автоматического выбора шаблонов или вызов данных из вышестоящих систем (опционально) - считыватели штрих-кода / QR-кода упрощают рабочие процедуры для оператора, а также предотвращают ошибки в работе.

LOGIN VIA NFC

The Qpix Control 2 software supports user login using an external NFC reader. Depending on the NFC tag/card, the customer's existing access cards can also be programmed in, for example.

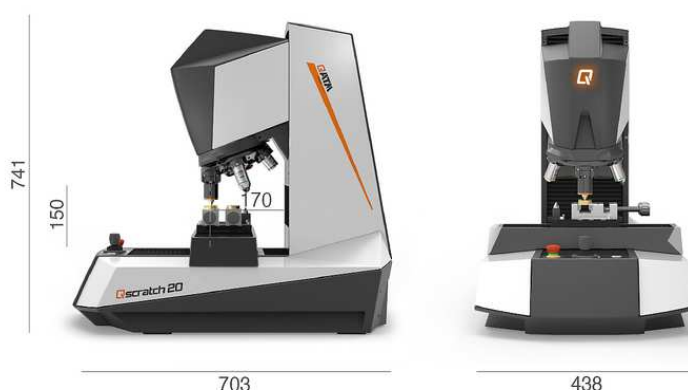

Диапазон испытательных нагрузок

Constant load

Progressive load

Z-ось

Позиции турели

Included basic equipment

Lens types

Обзорная камера

Cross slide

0.5* - 200 N

Single and multiple measurements

Single and multiple measurements

Z 150 mm (5.91")

8-fold motorized tool changer (max. 7 lenses)

Rockwell diamond; 5x, 20x, 50x Lenses

Standard (Achromat) for hardness testing and microscopy

5 Mpixel - CMOS color, USB 3.0; 52 x 39 mm (2.05 x 1.54")

Motorized X/Y axis

Размер предметного стола	150 x 120 mm (5.91 x 4.72")
Positioning accuracy	X/Y-axis +/- 0.2 µm
Диапазон перемещения	X 150, Y 150 mm (5.91 x 5.91")
Control elements	Emergency stop, Start button, Joystick X/Y/Z
Максимальный вес образца	50 kg (110 lbs)
Вес базовой машины	60 kg (132.3 lbs)
Optional	Hardness Upgrade and Qpix Inspect Module

* 0.5 N and 1 N > 1% accuracy

www.qatm.com/qscratch-20

ИНФОРМАЦИЯ ДЛЯ ЗАКАЗА

FULLY AUTOMATIC SINGLE MACRO/MICRO SCRATCH TESTER

BASIS MACHINE

QB0180101



Qscratch 20 GO

QB0180201

Qscratch 20

OPTIONS FOR QSCRATCH 20 GO

QA0117300

X Linear Slide 55 mm - Qscratch 20 GO

SCRATCH TEST UNITS AND TEST TOOLS

QA0103800

Scratch test module - Normal force - included

QER14000EA

Diamond Rockwell scratch tester 120° R0.2 - included

QER14000EU

Diamond Rockwell scratch tester 120° R0.2 - no cert

QSH14005EU

Scratch stylus carbide 0.5 mm, ISO 1818-1

QSH14010EU

Scratch stylus carbide 1.0 mm, ISO 1818-1

LENSES FOR SCRATCH TESTING AND QPIX INSPECT

QA0115900

Lens 10x/0.3 - FL 60 / WD 16

QSO220800

Exchange Lens 20x to 10x/0.25 FL 60 / WD 13

QA0116200

Lens 100x/0.90 - FL 60 / WD 1.0

HQ LENSES FOR SCRATCH TESTING AND QPIX INSPECT

QSO213300

Upgrade to HQ Lens 5x/0.15 - FL 60 / WD 23.5

QSO166400

HQ Lens 10x/0.30 - FL 60 / WD 17.5

QSO213400

Upgrade to HQ Lens 20x/0.45 - FL 60 / WD 4.5

QSO213500

Upgrade to HQ Lens 50x/0.80 - FL 60 / WD 1.0

QSO166700

HQ Lens 100x/0.90 - FL 60 / WD 1.0

WEAR PLATES FOR XY- SLIDE

QA0117800 Wear plate 150x120x8 mm - Qscratch 20 - included

HARDNESS TESTING OPTION

QA0117900 Vickers Hardness - Upgrade package "Qscratch 20"

QEV05000EA Test diamond Vickers - included (Knoop/Brinell/Rockwell also possible)



PC SYSTEM FOR MACHINE CONTROL

QA0057500 High-end PC system Full HD - Qpix Control2

QA0057600 2nd 24" Full HD UltraSharp monitor

QA0057501 High-end PC system Ultra HD (4K) - Qpix Control2

QA0057601 2nd 27" Ultra HD (4K) monitor

ДЕРЖАТЕЛЬ ОБРАЗЦА

QA0116300 Universal Precision fine vise

REFERENCE MATERIAL

QA0118000 QATM Scratch Reference Material - PACVD TiN 2–4 µm - included

SOFTWARE MODULES

QSO013600 Qpix - Barcode / QR code scanner software module

QA0089600 Barcode/QR-Code 1D/2D Scanner Kit - Bluetooth

QA0071700 El. Measurement caliper - Qpix Control2 compatible

QA0108200 NFC Login Package QpixControl2

QSO213600 RFID Tag MIFARE, key chain gray 10 pcs.

QPIX INSPECT - MICROSTRUCTURE ANALYSIS

QSO171100 Qpix Inspect - Phase analysis

QSO171200 Qpix Inspect - Layer thickness measurement

QSO171300 Qpix Inspect - Grain size determination

QSO194100 Qpix Inspect - Weld seam measurement

QPIX CONNECT - DATA CONNECTION

QSO031600	Qpix Control2 - PCI plug-in
QSO013500	Qpix Control2 - AMS plug-in
QSO030400	Qpix Control2 - Q-DAS plug-in

ПРИНАДЛЕЖНОСТИ

QSO165900	Dust protection cap - Qness 60 EVO, CHD Master+
A5810431	Vibration-Damped Lab Top
M7500006	Laboratory table 880x770x800 with stone plate
M7500009	Laboratory table 880x900x800 with stone plate
M7500007	Laboratory table 1320x770x800 with stone plate
M7500010	Laboratory table 1320x900x800 with stone plate
M7500008	Laboratory table 1760x770x800 with stone plate
M7500011	Laboratory table 1760x900x800 with stone plate
QA0048500	UPS - Uninterruptible Power Supply 900VA
QSO232600	Accessory box - Scratch tester
QVV0030	Documentation
QVV0001	Work certificate
QVV0002	Hardness comparison chart